

(19) World Intellectual Property Organization
International Bureau



(43) International Publication Date
9 January 2003 (09.01.2003)


PCT

(10) International Publication Number
WO 03/001972 A2

- (51) International Patent Classification⁷: **A61B**
- (21) International Application Number: PCT/US02/20158
- (22) International Filing Date: 27 June 2002 (27.06.2002)
- (25) Filing Language: English
- (26) Publication Language: English
- (30) Priority Data:
60/302,069 29 June 2001 (29.06.2001) US
60/372,579 12 April 2002 (12.04.2002) US
- (71) Applicant: **ETHICON, INC.** [US/US]; u.s. Route 22, Somerville, NJ 08876 (US).
- (72) Inventors: **TRACEY, MICHAEL, R.**; 403 Kacey Court, Branchburg, NJ 08876 (US). **KARL, Frederick, T.**; 27 Cooper Creek Circle, Newtown, CT 06470 (US).
- (74) Agents: **JOHNSON, Philip, S.** et al.; Johnson & Johnson, 1 Johnson & Johnson Plaza, NJ 08933 (US).
- (81) Designated States (*national*): AE, AG, AL, AM, AT, AU, AZ, BA, BB, BG, BR, BY, BZ, CA, CH, CN, CO, CR, CU, CZ, DE, DK, DM, DZ, EC, EE, ES, FI, GB, GD, GE, GH, GM, HR, HU, ID, IL, IN, IS, JP, KE, KG, KP, KR, KZ, LC, LK, LR, LS, LT, LU, LV, MA, MD, MG, MK, MN, MW, MX, MZ, NO, NZ, OM, PH, PL, PT, RO, RU, SD, SE, SG, SI, SK, SL, TJ, TM, TN, TR, TT, TZ, UA, UG, UZ, VN, YU, ZA, ZM, ZW.
- (84) Designated States (*regional*): ARIPO patent (GH, GM, KE, LS, MW, MZ, SD, SL, SZ, TZ, UG, ZM, ZW), Eurasian patent (AM, AZ, BY, KG, KZ, MD, RU, TJ, TM), European patent (AT, BE, CH, CY, DE, DK, ES, FI, FR, GB, GR, IE, IT, LU, MC, NL, PT, SE, TR), OAPI patent (BF, BJ, CF, CG, CI, CM, GA, GN, GQ, GW, ML, MR, NE, SN, TD, TG).

[Continued on next page]

(54) Title: A SYSTEM AND METHOD FOR ASSESSING URINARY FUNCTION

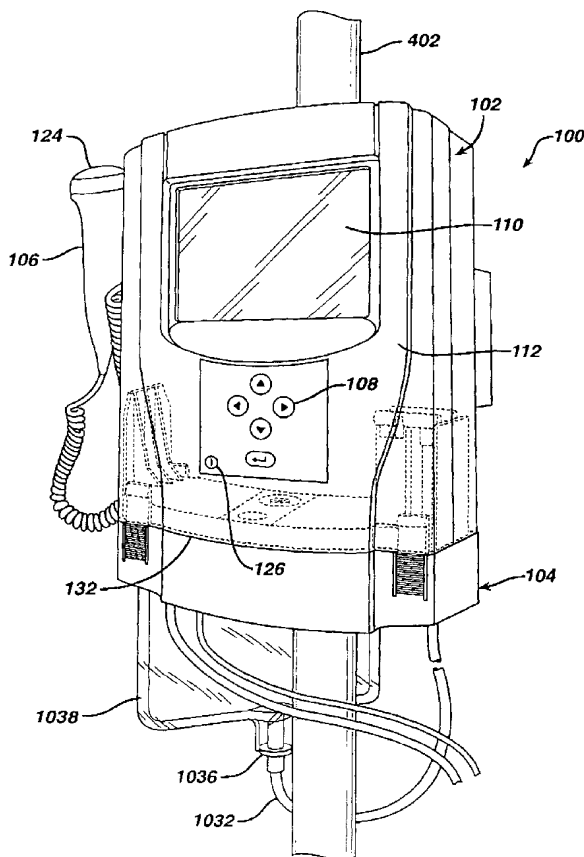

(57) Abstract: A vaginal speculum assembly is provided including an upper arm and a lower arm movably coupled to one another by a coupling member positioned at a posterior end. Each of the upper and lower arms are dimensioned for at least partial insertion into a patient's vagina, and when so inserted, are dimensioned to exert pressure on at least upper and lower portions of an inner surface of the patient's vaginal wall. The coupling member permits movement of the upper and lower arms relative to one another to adjust their positioning relative to the vaginal wall. The assembly further includes a locking device coupled to at least one of the upper and lower arms which, when activated, substantially prevents movement of the upper and lower arms relative to one another. A connector member extends outwardly from a posterior end of the upper arm. The connector member is configured to engage a pre-determined portion of an insert device assembly including an insert member to be inserted into the patient's urethral canal, to thereby maintain the insert device assembly in position relative to the vaginal speculum assembly.



WO 03/001972 A2



Published:

— without international search report and to be republished
upon receipt of that report

For two-letter codes and other abbreviations, refer to the "Guidance Notes on Codes and Abbreviations" appearing at the beginning of each regular issue of the PCT Gazette.

A SYSTEM AND METHOD FOR ASSESSING URINARY FUNCTION

Cross Reference to Related Applications

5

The present invention claims the benefit of earlier filed United States provisional patent applications, serial numbers 60/302,069, filed on June 29, 2001 and 60/372,579, filed on April 12, 2002, which are both incorporated by reference in their entirety herein.

10

Field of the Invention

The present invention relates generally to a system and a method for assessing urinary function. More particularly, the system and method is used for testing the integrity of the urinary system for diagnostic purposes and for use with therapies to correct urinary incontinence.

15

Background of the Invention

20

Women account for more than 11 million of incontinence cases. Moreover, a majority of women with incontinence suffer from stress urinary incontinence (SUI). Women with SUI involuntarily lose urine during normal daily activities and movements, such as laughing, coughing, sneezing and regular exercise.

SUI may be caused by a functional defect of the tissue or ligaments connecting the vaginal wall with the pelvic muscles and pubic bone. Common causes include repetitive straining of the pelvic muscles, childbirth, loss of pelvic muscle tone and estrogen loss. Such a defect results in an improperly functioning urethra. Unlike other types of incontinence, SUI is not a problem of the bladder.

25

Normally, the urethra, when properly supported by strong pelvic floor muscles and healthy connective tissue, maintains a tight seal to prevent involuntary loss of urine. When a woman suffers from the most common form of SUI, however, weakened muscle and pelvic tissues are unable to adequately support the urethra in its correct position. As a result, during normal movements when pressure is exerted on the

30

bladder from the diaphragm, the urethra cannot retain its seal, permitting urine to escape. Because SUI is both embarrassing and unpredictable, many women with SUI avoid an active lifestyle, shying away from social situations.

5 SUI is categorized into three types. Type I and Type II are directed to urethral hypermobility. Type III is directed to intrinsic sphincter deficiency (ISD). Diagnosis of ISD requires urodynamic evaluation. Urodynamic evaluation involves complex and invasive equipment and often requires referral to a specialist trained in urodynamic evaluation.

10 Existing diagnostic systems all require a catheter be passed trans-urethrally to measure pressure, such as Leak Point Pressure (LPP) –or Urethral Pressure Profile (UPP). An exemplary system is disclosed in publication (WO 0023127). Detection of LPP requires that a pressure sensor and catheter be passed trans-urethrally. The bladder is filled, and pressure is recorded. Fluid leakage from the urethral opening (meatus) corresponds to the maximum pressure the urethral sphincter can resist, or LPP. During
15 the UPP measurement procedure a pressure sensor tipped catheter is placed trans-urethral into the bladder and then withdrawn at a constant velocity. The pressure profile along the urethra, from bladder neck to meatus is recorded.

Other parameters may also be measured, such as abdominal pressure and urinary flow. A cystometrogram (CMG) is a pressure study that simultaneously
20 measures intra-abdominal, total bladder, and true detrusor pressures. Uroflometry measures urine flow rate visually, electronically, or via a disposable system. Video Urodynamic Systems also exist that simultaneously measure parameters, as described above, with radiographic visualization of the lower urinary -tract.

Existing urodynamic evaluation systems are complex, expensive, and require
25 extensive training. Furthermore, existing urodynamic systems often require at least 30 minutes to complete a test. This exceeds the time available for most standard physician office visits and results in referral to a specialist. No urodynamic system exists that can quickly and inexpensively record useful urodynamic measures, without passing a catheter or instrument trans-urethrally.

30 There remains a need for an improved system and method for assessing urinary function.

Summary of the Invention

A vaginal speculum assembly is provided including an upper arm and a lower arm movably coupled to one another by a coupling member positioned at a posterior end. Each of the upper and lower arms are dimensioned for at least partial insertion into a patient's vagina, and when so inserted, are dimensioned to exert pressure on at least upper and lower portions of an inner surface of the patient's vaginal wall. The coupling element permits movement of the upper and lower arms relative to one another to adjust their positioning relative to the vaginal wall. The assembly further includes a locking device coupled to at least one of the upper and lower arms which, when activated, substantially prevents movement of the upper and lower arms relative to one another. The assembly also includes a connector member extending outwardly from a posterior end of the upper arm. The connector member is configured to engage a predetermined portion of an insert device assembly including an insert member to be inserted into the patient's urethral canal to thereby maintain the insert device assembly in position relative to the vaginal speculum assembly. In one embodiment at least a portion of the connector member is comprised of a flexible material, and in yet another embodiment, the flexible portion can be manually manipulated to engage and grasp the insert device assembly.

According to yet another embodiment, when the upper and lower arms are inserted into the patient's vagina, the connector member is further configured to maintain the insert device assembly in position relative to the vaginal speculum assembly so that the insert member extends into the patient's urethral canal. In another embodiment, the insert member is a meatus plug member, and in an alternate embodiment, the insert member is a catheter.

A method is also provided for assessing urinary function including the steps of positioning a vaginal speculum assembly relative to a patient's vagina so that upper and lower arms of the vaginal speculum assembly support at least upper and lower portions of an inner wall of the patient's vagina, locking the upper and lower arms in position relative to one another, inserting an insert member of an insert device assembly at least partially into the patient's urethral canal, coupling the insert device assembly to the vaginal speculum assembly so that the insert member remains positioned within the

patient's urethral canal, and performing a urodynamic assessment utilizing the insert device assembly.

5 A vaginal speculum assembly is also provided including an upper arm and a lower arm dimensioned for at least partial insertion into a patient's vagina and movably coupled to one another at a posterior end. The upper and lower arms are movable relative to one another between a first collapsed position and a second open position, and, when inserted into the patient's vagina and in the open position, the upper and lower arms exert pressure on upper and lower portions of an inner surface of the patient's vaginal wall. The assembly further includes a locking device coupled to at least one of the upper and lower arms that, when activated, substantially prevents movement of the upper and lower arms relative to one another, and a connector member extending outwardly from a posterior end of the upper arm. The connector member is configured to engage a predetermined portion of an insert device assembly including an insert member to be inserted into the patient's urethral canal to thereby maintain the insert device assembly in position relative to the vaginal speculum assembly.

15 These and other features and advantages of the present invention will become apparent from the following more detailed description, when taken in conjunction with the accompanying drawings which illustrate, by way of example, the principles of the invention.

20

Brief Description Of The Drawings

Figure 1 is a perspective view of a one embodiment of a portable medical system according to the present invention;

25 Figure 2 is a front perspective view of a control device according to the present invention;

Figure 3 is a rear perspective view of the control device of Figure 2;

Figure 4 is a front elevational view of a control device in accordance with the present invention attached to a pole;

30 Figure 4a is an exploded perspective view of one embodiment of a pole attachment mechanism;

Figure 4b is a rear perspective view of the pole attachment mechanism of Figure 4a;

Figure 5 is an exploded perspective view illustrating interaction of a control device identification mechanism and module identification components;

5 Figure 5a is a schematic cross-sectional view taken across line 5a-5a of Figure 5 prior to engagement of the control device with the test module;

Figure 5b is a schematic cross-sectional view similar to Figure 5a showing engagement of the control device with the test module;

10 Figure 6 is a front perspective view of a module according to the present invention;

Figure 7 is a schematic illustration of one embodiment of control device electronics assembly;

Figures 8a-8i are flow diagrams illustrating operation of control device software and graphical user interface components;

15 Figure 9 is an alternate embodiment of a medical system according to the present disclosure;

Figure 10 is a schematic representation of a portable medical system including an SUI module;

20 Figure 10a is a partial cross-sectional view of one embodiment of a portable medical system including an SUI module;

Figure 11a is a side elevational view and partial cross-section of one embodiment of a hand actuator in an assembled configuration;

Figure 11b is a side elevational view and partial cross-section of the hand actuator of Figure 11a in an unassembled configuration;

25 Figure 11c is a side elevational view and partial cross-section of the hand actuator of Figure 11a in an operational mode;

Figure 11d is an alternative embodiment of a hand actuator according to the present invention;

30 Figure 12 is an enlarged perspective view of one embodiment of a meatus plug device;

Figure 13 is a schematic view of illustrating one embodiment of a urodynamic system in relation to a female urinary/reproductive system;

Figure 14 is a schematic view illustrating internal components of one embodiment of a system including a SCMG module;

Figures 15-16 are schematic views of the system of Figure 14 in relation to a female urinary/reproductive system;

5 Figure 17 is a schematic view illustrating one internal components of one embodiment of a system including a CCMG module;

Figures 18-19 are schematic views of the system of Figure 17 in relation to a female urinary/reproductive system;

Figure 20 is a flow diagram illustrating steps for using the system of Figure 10;

10 Figure 21 is a flow diagram illustrating steps for using the system of Figure 14;

Figure 22 is a flow diagram illustrating steps for using the system of Figure 17;

Figure 23 is a perspective view of one embodiment of an input pendant according to the present invention;

15 Figure 24 is a schematic view illustrating internal components of one embodiment of a system including a Uroflowmetry module;

Figure 25 is a schematic view illustrating use of the system of Figure 24;

Figure 26 is a perspective view of one embodiment of a vaginal speculum assembly in accordance the present invention;

20 Figure 27 is an exploded perspective view of the vaginal speculum assembly of Figure 26;

Figure 28 is a schematic view of one embodiment of a urodynamic system and speculum assembly in relation to the female urinary/reproductive system; and

Figure 29 is an exploded perspective view of a battery charger module that can be used in conjunction with the control device.

25

Detailed Description Of The Embodiments

30 Figures 1 through 26 illustrate generally various systems and methods for assessing urinary function and/or components of such systems and methods. Although the systems and methods disclosed herein are described in detail in relation to the female urinary system, it is to be understood that the present invention can readily be adapted for use in assessing male urinary function as well. Further, those skilled in the art will recognize that inventive principles, apparatus and methods disclosed herein may

also have application to assessing function in other areas, such as coronary function or pulmonary function. The present invention is to be limited only by the claims set forth herein.

Referring now to Figures 1 and 2, one embodiment of a portable medical system 100 is illustrated having particular application for assessing urinary function. The system 100 includes a control device 102 that controls operation of the system, at least one module 104 that can be removably coupled to the control device, at least one input device, such as the illustrated input pendant 106 and/or keypad 108, and at least one output device, such as the illustrated display screen 110. As will be described in more detail below, the control device 102 is designed to be removably coupled to any one of a plurality of testing modules 104 at any given time. As each module is uniquely suited to support a different type of diagnostic test or medical procedure, the resulting diagnostic system is not only readily portable, but is also extremely versatile in that the single control device, in conjunction with a plurality of small test modules, is capable of performing an array of diagnostic tests or other procedures. The system has particular application useful for assessing urinary function in that it provides a portable, modular system in contrast to the non-portable, expensive, and cumbersome equipment that is currently used for assessing urinary function. In addition, as will also be described in greater detail below, the present invention can perform tests quicker, and in a manner that is less uncomfortable and less invasive for a patient.

The control device 102 includes a housing 112 for housing various components, including one or more batteries 114, an electronics assembly 116, a pump device 118 including a motor, and various other circuitry. Batteries supply power to the control device 102, and are contained within a battery compartment 120 that is accessible by removing the battery cover 122 that forms part of the housing 112. The control device further includes an input keypad 108 for allowing a user to input data (such as patient name or other identifier, numeric identifiers, patient history, date etc.) and an input pendant 106 including one or more switches 124 that allow user input of additional information (i.e., event input based on patient feedback), and an activation switch 126 for turning the device on and off. The pump device 118 and at least one pressure transducer 128 are also contained within the housing. The pump device is electrically coupled to the battery and the electronics assembly, and the pressure transducer is electrically coupled to the electronics assembly. The control device 102 may also

include a pole mounting mechanism 400 for mounting the control device on a pole such as the pole of an IV solution caddy 402 including a hook 404. One embodiment of a pole mounting mechanism is illustrated in Figures 4a and 4b. The device may also include an interface 130 including appropriate electrical pinouts to enable the control
5 device to communicate for purposes of battery recharging or printing of patient test data.

As indicated above, any one of a plurality of modules 104, such as diagnostic test modules, can be removably coupled to the control device 102, and the control device is designed to uniquely identify the attached module, and perform routines
10 specific to that module. Thus, the control device includes a module detection mechanism 500 capable of identifying the attached module that is electrically coupled to the electronics assembly (see Fig. 5). This module detection mechanism includes one or more identification probes 502 that project from the interface side 132 of the control device and are electrically coupled to the electronics assembly. The modules
15 104 may include one or more apertures in the module housing 506 that are designed to receive therein the identification probes when the module is removably coupled to the control device. When so coupled, the identification probes will bridge one or more module identification elements or components 504, such as resistors, capacitors, fuses or other suitable electronic components, present within the module. The identification
20 probes are electrically coupled to the electronics assembly 116 (described more fully below), which determines a value, such as resistance, associated with the module identification element(s) that they bridge. Each module is designed to have a value so that identification of this value by the electronics control assembly enables the control device to uniquely identify the attached module. In a preferred embodiment, the
25 control device may include one or more sets of identification probes 502 at different locations, and different modules have a module identification components 504 at different locations. The location, as detected by the control device, identifies the attached module. In yet another embodiment, the module identification component(s) may be coupled to an exterior side of the module housing so that apertures in the
30 module housing are not required.

The module further includes at least one coupling element 600 for removably coupling the module to the control unit (see Figure 6). In the illustrated embodiment, the module includes four coupling elements placed toward the ends of each of the front

and rear faces 602, 604 of the test module. Each coupling element contains a tab element 606 that engages a corresponding ridge 607 (best seen in Fig. 5) on an interior surface of the control device when the module is removably coupled to the control device. To couple the module to the control unit, the coupling elements are depressed slightly in the direction indicated by the arrow in Figure 6. The module is then aligned with the control device as shown in Figure 1, and the coupling elements released to allow engagement with the corresponding ridges described above. The module can subsequently be removed from the control unit by once again depressing the coupling elements and removing the module from the control device.

Finally, the module housing 506 includes first 608 and possibly second 610 ports therein as shown in Figure 6. Each of the first and second ports are configured so as to define a recess capable of receiving a control device pressure sensor, such as a pressure transducer, therein when the module is coupled to the control device. For example, a first control device pressure transducer 128 is received within the first port recess 608 and comes in physical contact with a pressure interface 1024 (see Fig. 10) so that pressure changes at the pressure interface can be transmitted to and detected by pressure transducer 128 and converted to electrical signals that are sent to the electronics assembly for interpretation. Similarly, the second port 610 also defines a recess capable of receiving therein a second control device pressure transducer 1030. The first and second ports are further configured to form an airtight seal with the control device when coupled thereto, preferably by incorporating sealing elements such as gaskets or the like. Individual modules and their operation in conjunction with the control device will be described in greater detail below.

As indicated above, contained within the housing 112 of the control device 102 is an electronics assembly 116 (see Figure 7) that is designed to control operation of the pump device 118, to acquire and format data from the pressure transducer(s), to drive a display 110 and/or other output device, and to accept and interpret input data, such as from switches 108, 126, and/or 124. The electronics assembly 116 consists of an integrated circuit board 702, hardware interfaces to the pump device 708, pressure transducer 706, 707, display 709 and switches 703, 704 and 705; and a microprocessor 710. The microprocessor 710 serves as the main controller for the diagnostic system and is supported by the custom integrated circuit 702 and powered by the batteries. Also included are interface connection elements including an electronic module

identification connection 712 to the electronic detection mechanism 500, and electronic connections 714 that enable downloading of data to a printer or other external device.

The microprocessor 710 is programmed with a custom program file. In the illustrated embodiment, this software has multiple functions. First is the acquisition of
5 input from the operator. This input data is captured from the input keypad 108, and/or switches 124, 126, pressure transducer(s) or other input device, depending upon which test module is in use. The software also controls operation of the pump device 118. Input data is interpreted and appropriate signals are sent to the pump device motor via the integrated circuit board 702. Yet another function is to acquire and condition data
10 from the pressure transducer(s). This data is then sent in the appropriate format to the display 110, along with applicable pump device data in the form of volume or time information. Finally, as indicated above, the software receives input from the module detection mechanism 500 and interprets this input to determine which test module is coupled to the control device.

15 Figure 8a-8i are flow diagrams illustrating operation of the diagnostic system software and features of the system graphical use interface for a preferred embodiment of the invention. When the system is powered on, the user is first presented with a welcome screen. While this screen is being displayed the system is undergoing a self-test routine 802 to test the integrity of system hardware and software components.
20 Upon completion of this routine, the user is provided with information relating to the amount of available system memory 804. Following the pressing of any key 806 on input device 108 by the user, the system identifies the attached module 808 as described above, and following such identification, the processor executes a software subroutine specific to the identified module. For each software subroutine, however, a main menu
25 is displayed next, such as that indicated by reference numeral 810. In the illustrated embodiment, the main menu includes six possible selections. "Utilities" enables the user to access various system features, such as setting the date, time etc, or adjusting the brightness or contrast of the screen; "Quit" terminates the session; "Patients" enables the user to access any previously stored data relating to other patients and tests already
30 performed; "Prime" initiates the pump priming process; "Patient ID" enables the user to enter a patient identification number; and "Test" initiates a software subroutine specific to the attached module to carry out the desired test procedure. In the presently described embodiment, the software and user interface associated with the "Prime," "Utilities,"

“Quit,” and “Patient ID” selections are substantially the same for each software subroutine. The “Test” and “Patients” selections, however, are different for each test module. Each of these selections will be described in greater detail below.

As is illustrated in Fig. 8a, the first time the main menu is displayed both “Test” and “Prime” appear in a different color or shade from the other options, indicating that they are not currently available. This is to ensure that patient identification information is entered before proceeding with any priming or testing procedures. The user may select the “Patient ID” option by scrolling using the appropriate arrows on the input keypad 108. Following this selection the Patient ID screens appears 820 (Fig. 8b). In the illustrated embodiment, the patient ID consists of a nine digit integer. To enter the patient ID, the user scrolls to a selected blank using the left and right arrows and/or left and right arrows on the input keypad 108 (824) to select desired numbers. Once the desired number is selected, the user presses ENTER; the selected number will then appear in the rightmost blank. Subsequent numbers are selected as described above, and will appear in the rightmost blank while previously selected numbers move to the left. This process is completed until all blanks are filled in. In one embodiment, there is a default value for each blank, such as 0, and the user may proceed with testing by accepting the default patient ID number consisting of all 0’s. Once complete patient identification information is entered, the user selects the “Main Menu” option 832, which returns to the main menu screen. At this point, however, the “Prime” option become available 834 (and “Patient ID” is no longer available).

Before performing any test that requires fluid to be infused into the patient, priming operations must be performed to ensure that the fluid infusion lines (tubing) are filled with fluid and not air. Referring now to Fig. 8c, the user selects the “Prime” option 840 by using the arrow keys to select the option, and then pressing the enter key. The Prime screen then appears. According to one embodiment, the Prime screen includes two options as indicated at 842: “Prime” or “Main Menu.” In another embodiment, the Prime screen is particular to each module, and may present only one option to initiate priming. Selecting the Prime option causes the pump to start and run for a predetermined amount of time, such as 20 seconds, and then automatically shuts off. The user is then presented with a screen 846 at which the user can accept the prime as complete (MAIN), or choose to reprime (PRIME). When priming is accepted as complete, the main menu once again appears, this time with “Test” as an option 848. In

another embodiment, priming operations may be specifically tailored for different test modules. For example, as will be described in more detail below, the SUI test modules includes a hand actuator including an activation button 1118 or 1128. The system may be designed so that following display of the Prime screen, pump priming operations can
5 be initiated by depressing the activation button.

With priming complete, testing can begin. As indicated above, testing procedures depend on the attached test module, and accordingly, the software and graphical user interfaces relating to each test module will be discussed in greater detail below in conjunction with the detailed description of each test module.

10 In an alternative embodiment of the invention illustrated in Figure 9, the control device 102 is electrically coupled to a laptop/standard computer 900, and the microprocessor and associated software reside in the computer.

As indicated above, the diagnostic system described herein has particular application to urodynamics in that it enables clinicians to diagnose a plurality of urinary
15 incontinence problems when used with specifically designed testing modules (to be discussed hereinafter). As a miniaturized urodynamic tool, the control device 102 in conjunction with modules 104 can measure urethral resistance pressure (URP), voiding flow (Uroflometry), and bladder dysfunction (Cystometrogram (CMG)). As will be described further below, URP is a new and unique approach to urodynamic
20 measurement of stress incontinence that is less invasive for a patient, and faster than currently known and used diagnostic tests. Uroflometry is the study of micturation over time. CMG is the study of bladder or detrusor instability. A major advantage of the diagnostic system disclosed herein is that it can achieve all of the diagnostic tests described above with a portable unit that can be used in any office exam room,
25 removing the need for the reservation or scheduling of a specialized urodynamic room, and the need for the complex equipment currently required for such tests. The urodynamic system is easy to use and does not require advance training. Use of the disclosed system makes testing more comfortable for patients by enabling faster set up, shorter test time, and less invasive procedures.

30 In actual use, different modules can be removably coupled to the control device 102 to conduct these different urodynamic tests. Each module performs a different and distinct test. These modules include, but are not limited to, a stress urinary incontinence (SUI) module for measurement of urethral resistance pressure (URP); a

simple CMG module for measurement of bladder instability; a complex CMG module for measurement of bladder instability; and a uroflometry module for the study of micturation over time. Modules may be suitably adapted to either male or female incontinence diagnosis.

5 Before proceeding with a discussion of individual test modules, to assist the reader a brief overview of the female urinary system will be described with reference to Figure 13. The female urinary system 1300 includes an elongated urethral canal 1302 having a urethral meatus (entrance) 1304 and having a substantially circular-shaped urethral sphincter muscle 1306 attached thereto, and a bladder cavity 1308 surrounded
10 by a detrusor muscle 1310. The detrusor muscle 1310 also surrounds and supports the urethral canal 1302. The bladder cavity 1308 is in close proximity to the abdominal wall 1312, the pubis bone 1314, the pelvic floor 1316 (levator ani muscle), the vaginal canal 1318, the clitoris 1320, the uterus 1322 and the anal sphincter muscle 1324.

Individual testing modules will now be described in detail.

15

STRESS URINARY INCONTINENCE MODULE

Figures 10-13 illustrate one embodiment of a stress urinary incontinence testing module (SUI) 1000 for diagnosing the involuntary loss of urine during physical
20 activities such as coughing, sneezing, laughing or lifting. The SUI testing module 1000 includes a SUI module housing 1002 that can be removably coupled with the control device 102 as described above. The module housing may be in the form of a plastic disposable cartridge. Within the module housing is a tubing assembly 1004 including a fluid inlet 1006, a fluid outlet 1008, and a first fluid conduit 1010 extending
25 therebetween. Tubing loop 1012 forms part of the tubing assembly and is positioned so that, when the SUI testing module is coupled to the control unit, the stator 1014 of the pump device 118 in the control unit 102 cooperates physically with the tubing loop 1012 so that the pump device operates as a peristaltic pump to pump fluid through the first fluid conduit 1010. To assist in this regard, a tubing guide 599 aids in positioning
30 a portion of the tubing assembly so that it will properly and effectively engage the peristaltic pump. According to the illustrated embodiment, tubing guide 599 has a substantially U-shaped configuration, however, many other configurations are suitable, as the principles of operation of peristaltic pumps are well known in the art. Tubing

member 1050 also forms part of the first fluid conduit. The module housing 1002 also includes a pressure chamber 1016 for dampening pressure fluctuations that may be caused by operation of the pump device. The pressure chamber 1016 is in fluid communication with the first fluid conduit 1010 via valve openings 1018a-c of three-way valve member 1020. The pressure chamber is filled primarily with air, but varying amounts of fluid may also be present. Positioned at a distal end of pressure chamber 1016 is a filter component 1022 designed to isolate fluid from electronic elements of the system 100. In this regard, filter 1022 may be a hydrophobic filter that allows air to pass into pressure interface 1024, but not liquid. When the testing module is coupled to the control device 102, pressure interface 1024 is in physical contact with pressure transducer 128 of the control device so that pressure fluctuations within the pressure chamber 1016 and pressure interface 1024 can be transmitted to and sensed by the pressure transducer, and subsequently transmitted to the electronics assembly as indicated above. In this manner, the control device measures pressure within the first fluid conduit of the tubing assembly of the SUI testing module, which substantially corresponds to the pressure within the urethral canal as described more fully below.

The SUI testing module 1000 tubing assembly also includes a second tubing member 1025 having a channel therethrough forming a second fluid conduit between a proximal end 1026 and a distal end 1028.

Referring now to Figs. 11a-c, the SUI testing module may also include a hand actuator 1100 having and insert device such as a meatus plug device 1102 attached thereto. The meatus plug device 1102 (see Figure 12) includes an attachment member 1104 at a proximal end 1106 coupled to a plug or insert element or member 1108 at a distal end 1110, and a channel 1112 extending therethrough allowing fluid flowing through the first fluid conduit to flow through the meatus plug device. The distal end 1114 of the plug element may also include one or more transversely aligned apertures or openings 1116 therein approximately equally spaced apart from one another around the exterior surface of the distal end. As the outer diameter of the distal end at the location of the apertures is less than the diameter of the inner wall of the urethral canal at that location (described more fully below), one or more of the apertures 1116 can be used for assurance of fluid flow into the urethra during actual operation.

In one embodiment, the hand actuator further includes a hand-sized housing or casing 1102 including therein an initiator element 1118 (Figs. 11a-c) that is in fluid

communication with tubing member 1025. Preferably, initiator element is an air bladder 1097 coupled to a distal end 1028 of the tubing member 1025. The proximal end 1026 of tubing member 1025 coupled to a pressure interface 1026a that is positioned so that, when the SUI testing module is coupled to the control device, pressure within tubing member 1025 can be sensed by pressure transducer 1030. As a closed system, pressure on the activation button 1118 can be sensed at the pressure interface 1026a by pressure transducer 1030, and interpreted by control device 102 as a signal to initiate and/or deactivate the test.

The hand actuator 1100 further includes a fluid conduit 1050 extending between an outlet 1195 and an inlet 1194 that is coupled to (integrally or otherwise) an external tubing conduit leading to a fluid source, such as the first fluid conduit 1010 of the SUI test module. Alternatively, the hand actuator may be designed to include therein the fluid source. The fluid outlet 1195 is in fluid communication with the insert member channel of the meatus plug device. An activation device 1127 including a trigger 1128 extends through an opening 1118a to an exterior of the casing. The activation device 1127 is movable between a first rest position (shown) and a second activated position. In the first position spring 1130 exerts force on coupling member 1132, causing it to pivot relative to pivot element 953 and pinch the distal ends of at least tubing member 1050 to prevent fluid flow therethrough. When in the second position, movement of the trigger causes the coupling member 1132 to pivot to a point at which it no longer pinches tubing member 1050. Further, trigger 1128 may also compresses air bladder 1097 to initiate testing as described above in connection with initiator element.

The plug element 1108 is configured so that, when inserted into the urethral meatus of a patient (see Figure 13), it will substantially block or prevent fluid flow out of the urethra, as well as into the urethra other than through the meatus plug device channel 1112. Further, when inserted, the plug element is positioned distal of the urethral sphincter 1306 (toward the outside of the body) as shown in Fig. 13. In the embodiment shown in Fig. 12, the distal end or distal portion 1114 of the plug element is substantially conical in shape, and decreases in diameter toward its distal end 1114. A proximal portion 1199 is configured to engage the inner wall of the urethral canal to substantially prevent fluid flow therebetween. Other shapes, however, are also possible so long as fluid flow into or out of the urethral is substantially blocked (other than through the meatus plug device channel) and the plug element remains located distal of

the urethral sphincter. The meatus plug device 1102 is made of a biocompatible material, such as stainless steel or polypropylene. The meatus plug device may be disposable, but may also be made of a sterilizable material so that it can be reused.

5 The first fluid conduit 1010 of the tubing assembly also includes an elongated single lumen tubing member 1032 having a first end 1006 and a second end 1034 and a fluid channel extending therethrough. A spike device 1036 is coupled to the first end 1006 of the single lumen tubing member for attachment to a fluid bag 1038 (having a fluid 1010 therein) in a manner well known in the art. As described above, the meatus plug device and first fluid conduit are coupled to one another such that fluid from the
10 fluid source traveling through the first fluid conduit may pass through the insert member (via the channel therein) and into the urethral canal distal of the urethral sphincter. Further, as the first pressure interface 1024 is in fluid communication with the first fluid conduit and ultimately the urethral canal, pressure at the pressure interface substantially corresponds to the pressure within the urethral canal distal of the urethral
15 sphincter.

Use of the system 100 including a SUI testing module 1000 is as follows. First, the SUI testing module is removably coupled to the control device 102 in the manner described above. The physical coupling causes the identification probes 502 of the control unit to engage the module identification element(s) 504 of the SUI testing
20 module, enabling the control device to identify the SUI testing module. The physical coupling also brings pressure interface 1024 in physical contact with pressure transducer 128 as described above so that pressure changes at the pressure interface can be detected by the pressure transducer and transmitted to the electronics assembly for interpretation. The pressure interface 1026a at the proximal end of tubing member
25 1025 similarly comes in contact with pressure transducer 1030 so that pressure within tubing member 1025 can also be detected. Finally, the tubing loop 1012 is brought into physical contact with the pump device 118 so that the pump device can drive fluid through the first fluid conduit by peristaltic motion, as described above.

As shown in Fig. 20, once the SUI testing module 1000 is coupled to the control
30 device 102 (2010), the operator enters appropriate input data into the keypad 108 or other input device (2015) for the SUI test (described in more detail below). This data is received and interpreted by the microprocessor 710 and applicable information is sent by the microprocessor to the display 110. Priming operations are then performed

(2020) to ensure that the first fluid conduit 1010 contains fluid. At this point, the microprocessor is ready to start the test routine.

The meatus plug 1102 is inserted into the meatus of the urethra (2025) and the test is started (2030) by pressing the activation button as described above. This in turn
5 sends instructions to the pump device via the integrated circuit. The pump device then pumps fluid 1040 through the first fluid conduit 1010 and meatus plug device channel 1112 and into the urethral canal distal of the urethral sphincter (2035). As fluid pressure builds in the urethral canal 1302, pressure in the pressure chamber 1016 also builds. This pressure is transmitted through the filter component 1022 and pressure
10 interface 1024 to the pressure transducer 128, which receives the pressure data and transcribes it into an electrical signal. The electrical signal from the pressure transducer is sent to the microprocessor 710 via the integrated circuit 702 where it is acquired and conditioned. The information is then sent to the display 110 via the integrated circuit. The microprocessor ends the test after a specified amount of time, or upon receipt of
15 input from the user by sending an “off” signal to the pump motor drive. Once the test has been completed, the operator disengages the activation button 1118 (step 2040) and removes the meatus plug element from the meatus 1304 (2045).

Referring once again to Figures 8a-i, and in particular Fig. 8d, when the “Test” option is selected the SUI test can be performed. The SUI Test screen appears 860, and
20 the user initiates the test by depressing the trigger 1128 or movable shell 1126 (862) to allow fluid flow into the urethral canal as described above. The motor is then activated and the pump device pumps fluid into the urethral canal for a predetermined period of time, preferably 15 to 20 seconds. During this time a graph (see 860) is continuously displayed illustrating measured pressure on the vertical axis (preferably in cm of water)
25 versus time on the horizontal axis. As fluid is pumped into the urethral canal, pressure within the urethral canal distal of the sphincter continues to increase until that point in time at which the urethral sphincter yields (open) under the force of the pressure within the urethral canal. At that point the pressure curve becomes substantially flat, as illustrated in Fig. 8d, since the sphincter is open and fluid is filling the bladder. The
30 value of the flat portion of the curve is considered the “urethral resistance pressure (URP),” and can be obtained from the displayed graph. On completion of the test (after expiration of the predetermined time period the pump device stops), the graph remains, and the user is preferably provided with an option to adjust the software generated URP

value (860a) before saving the test results. To adjust the URP value, the user uses the up and down arrows to manipulate a horizontal line which indicates the URP value that appears on the screen (870). When the ghost line is at the desired value, the user presses enter (872).

5 Once the final URP value is displayed, a Save/Delete screen 874 is overlaid on the screen. If the user selects the "Save" option, the test results are saved in memory. If the user selects "Delete" from the Save/Delete screen 874, the user is then presented with the Save Test screen 876. If "Delete" is chosen the test is deleted, but if "Cancel" is selected, the user is returned to the Save/Delete screen.

10 According to one embodiment, test results for up to three out of six possible tests may be stored. Once three tests have been stored or six tests have been run, whichever comes first, the control unit 102 will disable the module identification component 504 via the identification probes 502. After testing is complete, the user may return to the main menu by selecting the "Menu" option from the Test Complete
15 screen.

 One option available from the Main Menu, as stated above, is "Patients," which allows the user to access patient and test data previously stored. According to one embodiment illustrated in Fig. 8h, when "Patients" is selected from the Main Menu, a Patients Screen 891 appears. On this screen, options for each patient and test for which
20 data has been stored 892 are presented and selection of one of these options causes a Patient Test Menu 893 to be displayed (Fig. 8i). Selecting "Delete" 896 will present the user with the option to delete the stored data for that patient/test, and selecting "Print" 895 will enable the user to print the stored data. The Print option will only be available (will not be greyed out) when the control device is coupled to a cradle, or
25 otherwise appropriately coupled to a printer. Selecting "View Test" will cause a Patients Test screen 898 or 899 to appear depending on whether stored data is a CMG (898) or a SUI (899) data set. The Patients Test screen may vary depending on the test module that is attached. For example, for the SUI stored data, the Patients Test screen is the screen illustrated by 899, whereas for the CMG data (discussed below), the
30 Patients Test screen is the screen illustrated by 898. The Patients Test screens provide the user with the option to view data relevant to the particular form of test performed.

 As indicated above, the results obtained from the SUI test is the urethral resistance pressure (URP), which is the back-pressure necessary to force open the

urethral sphincter muscle 1306 from the reverse or opposite direction from which fluid normally flows. A major advantage of the SUI testing module 1000 is that the insert or plug element 1108 of the meatus plug device 1102 only enters the external urethral canal (meatus) and does not cause any discomfort associated with passing a catheter through the internal urethral sphincter. Thus, the diagnostic system disclosed herein having a SUI module 1000 is less invasive and more comfortable for patients. Further, the testing procedure for the SUI module 1000 is easy to implement, quick to perform, and does not require advance training by the clinician and/or physician.

10 SIMPLE CYSTOMETROGRAM (CMG)

The diagnostic system disclosed herein can also be used to perform both simple and complex cystometrograms. Figures 14-19 show both simple (SCMG) and complex cystometry (CCMG) systems for the testing of bladder function in which pressure and volume of fluid in the bladder cavity 1308 is measured during filling, storage and voiding. Urologists typically measure the static pressure relationship in the bladder of patients, this being termed as a cystometrogram (CMG), in order to determine the capacitance of the bladder as a function of pressure and volume.

Referring now to Figure 14, the SCMG testing module 1400 includes a module housing 1020b that can be removably coupled to the control device 102 in the manner described above. The module housing 1020b may be in the form of a plastic disposable cartridge. The SCMG testing module contains many elements that are similar to those described above in connection with the SUI testing module, and thus like numerals will be used for these elements. Contained within the module housing is tubing assembly 1004b including a first fluid conduit 1402 between fluid inlet 1404 and fluid outlet 1406. The tubing assembly also includes a second conduit 1408 between a distal end 1410 and a proximal end 1412. Coupled to the proximal end is a filter 1022b and pressure interface 1412 that contacts pressure transducer 128 to convey pressure information thereto when the SCMG testing module is coupled to the control device. Compliant tubing loop 1012 similarly forms part of the first fluid conduit, and couples with the pump device 118 in the same manner as described above in connection with the SUI module. The distal ends 1406, 1410 of the first and second conduits are each coupled to respective proximal ends 1414, 1416 of first and second tubing elements

1418, 1420 of a dual lumen catheter 1422 so that the first and second conduits 1402, 1408 between the proximal 1414, 1416, and distal 1460, 1462 end of the dual lumen catheter are in fluid communication with channels in the first and second tubing elements 1418, 1420 of the dual lumen catheter 1422. This attachment may be
5 accomplished by an adhesive bond, a solvent bond, an ultrasonic weld, or any other suitable type of attachment that creates a fluid tight seal. In another embodiment, the dual lumen catheter is an inflatable balloon catheter such as a Foley-type catheter, that includes a pressure sensor 1424 positioned at the tip of the catheter (see Figure 16). Any other suitable catheter may also be used, such as fiber optic or air charged
10 catheters. The pressure sensor may be a micro tip transducer, an air charged sensor, a fluid charged sensor, a fiber optic sensor or any other pressure measuring sensor.

Use of the diagnostic system to perform a SCMG will now be described in detail with reference to Figures 15, 16 and 21. First, the SCMG testing module is coupled to the control device in the manner described above (2110). The physical
15 connection causes the identification probes 502 of the control unit to engage the module identification element(s) 504 of the SCMG testing module, enabling the control device to identify the SCMG testing module in the manner described above. The physical coupling also brings the pressure interface 1024b in contact with the pressure transducer 128 so that pressure changes in the second fluid conduit can be detected by
20 the pressure transducer. This coupling also causes the tubing loop 1012 to engage the pump device so that the pump can drive fluid through the tubing loop by peristaltic motion, as is also described above.

Once the SCMG testing module 1400 is coupled to the control device 102, the operator enters input data appropriate for the SCMG test (2115). This data is received
25 and interpreted by the microprocessor 710 and applicable information is sent by the microprocessor to the display 110. Priming operations are then performed (2120). At this point, the microprocessor is ready to start the test routine.

The dual lumen catheter 1422 is then inserted into the bladder 1308 (2125) via the urethra 1304 and the test is started by pressing the input pendant switches 124
30 (2130). The microprocessor 710 receives the signal from the input pendant switches. Instructions are then sent to the pump device 118 via the integrated circuit 702. The pump device then pumps fluid through the first fluid conduit 1402 and tubing element 1418 into the bladder (2135). As fluid volume builds in the bladder, pressure in the

bladder also builds. This pressure is transmitted through tubing member 1420 and the second conduit 1408, filter component 1022b, and pressure interface 1024b. The pressure transducer 128 receives the pressure data and transcribes it into an electrical signal. The electrical signal from the pressure transducer 128 is sent to the
5 microprocessor 710 via the integrated circuit board 702 where it is acquired and conditioned. During the course of a typical SCMG test, the patient provides event input, such as feeling the need to void and/or the intensity of that feeling, which is input to the control device via input pendant switches 124, as will be described more fully below. The microprocessor ends the test (2140) after a specified amount of time, or
10 upon receipt of an "off" signal from input pendant switch 124. Once the test has been completed, the operator removes the catheter 1422 from the bladder (2145). Following the test the software then exits the SCMG test subroutine, and the data storage routine is run to store and/or display results of the test.

Referring again to Figures 8a-i, and in particular Figure 8e, when the
15 "Test" option is selected the SCMG test can be performed. The SCMG Test screen appears 870a, and the user initiates the test by depressing input pendant switch 124 (see Fig. 10a) quickly. The pump device is then activated and pumping begins 872a. In a preferred embodiment, fluid is infused into the patient's bladder at a rate of approximately 1 ml/sec. As such, this test may be approximately 16 minutes in
20 duration, as opposed to approximately 15-20 seconds that may be required for the SUI test.

As the bladder is filling, the patient communicates the point in time at which he/she feels the initial sensation of needing to void, and the user presses the input pendant switch 124 to mark this point in time 873a. The fluid infusion continues, and
25 the user then marks the point in time at which the patient feels the urge to void 874a, and the point at which the patient feels an extreme, almost unbearable urge to void 875a, or has voided. Upon this third marking, the fluid infusion ceases and the test is completed 876a. During fluid infusion and after the test is complete, a graph is displayed illustrating pressure versus volume infused. After completion of the test a
30 Save/Delete overlay 877 appears. Selecting "Save" and pressing enter saves the test data. Selecting "Delete" causes a Save/Delete screen 878 overlay to appear. Selecting "Delete" from this screen deletes the data, where as selecting "Cancel" from this screen returns to the Save/Delete overlay.

At any point between initiating pumping and completing the SCMG test, the user may pause the test by depressing and holding, or pressing firmly on the input pendant switch 880, which causes the pump device to stop pumping fluid into the patient's bladder, and a Pause screen 881 (Fig. 8f) to appear on the display. Selecting "Quit" causes a End Test screen 885 to appear, and if "OK" is selected the test is stopped 886. If "Cancel" is selected the Pause screen reappears. If "Resume" 883 is selected from the Pause screen 881, the SCMG test resumes where it left off (pumping begins again). If, however, "LPP" 882 is selected from the Pause screen 881, assessment of the patient's leak point pressure (LPP) begins. No pumping of fluid occurs during this test. First, a LPP screen 887 appears and a blank graph is displayed. Pressure in centimeters of water is plotted on the vertical axis versus time on the horizontal axis. The patient then proceeds to exert pressure on the bladder as if attempting to void 888. The user marks the point at which a leak occurs 889, and the test is automatically completed after three minutes or three leaks, upon which the user is returned to the Pause screen 881. LPP results may then be stored or deleted, the CMG test may be resumed, or the test can be terminated altogether.

COMPLEX CYSTOMETROGRAM

20

In reference to Figures 17-19, the complex CMG (CCMG) testing module 1700 is similar to the SCMG testing module, but the tubing assembly also includes an additional single lumen tubing member 1702 having a proximal end 1704 and a distal end 1706 and a third conduit extending therethrough. The proximal end 1704 of the single lumen tubing member is coupled to another filter component 1022c and pressure interface 1024c. Pressure interface 1024c contacts pressure transducer 1030 when the CCMG testing module is coupled to the control device, enabling pressure transducer 1030 to sense pressure within the third fluid conduit.

Use of the diagnostic system to perform a CCMG will now be described in detail with reference to Figures 18, 19 and 22. First, the CCMG module is coupled to the control device (2210). The physical connection causes the identification probes 502 of the control device 102 to engage the identification elements 504 of the CCMG testing module, enabling the control device to identify the CCMG testing module. The

physical coupling also brings pressure interfaces 1024b, 1024c in contact with the pressure transducers 128, 1030 so that pressure changes in the second and third conduits can be detected by the pressure transducers. This coupling also causes the tubing loop 1012 to engage the pump device 118 so that the pump can drive fluid
5 through the tubing in the CCMG module.

Once the CCMG testing module 1700 is coupled to the control device 102, the operator enters input data appropriate for the CCMG test (2215). This data is received and interpreted by the microprocessor 710 and applicable information is sent by the microprocessor to the display 110. Priming operations are then performed (2220
10 The dual lumen catheter 1422 is inserted into the bladder via the urethra 1302 (2225). The single lumen catheter 1702 is inserted into either the vagina or the rectum (2230) and the test is started (2235) by pressing the input pendant switches 124. The microprocessor 710 receives the signal from the input pendant switches. This in turn sends instructions to the pump device 118 via the integrated circuit 702, and the
15 pump device pumps fluid through the first tubing conduit 1042 and tubing element 1418 into the bladder (2240). As fluid volume builds in the bladder, pressure in the bladder also builds. This pressure is transmitted through pressure interface 1024b to pressure transducer 128. Similarly, abdominal pressure is transmitted through pressure interface 1024c to pressure transducer 1030. The pressure transducers receive the
20 pressure data and transcribe it into electrical signals. The electrical signals are sent to the microprocessor 710 via the integrated circuit board 702 where it is acquired and conditioned. The microprocessor ends the test after a specified amount of time or upon receipt of an "off" signal from input pendant switches 124 (2245). Once the test has been completed, the operator disengages the input pendant switches and removes the
25 catheters 1422 and 1702 from the bladder (2250). The stored information is then available for review on the display screen, or by a print out through a charging cradle (printer assembly), or downloaded to a PC via a software interface in the charging cradle.

Referring again to Figures 8a-i, the CCMG module software subroutine
30 and graphical user interface is substantially as described in connection with the SCMG module. The system subtracts the abdominal pressure from the bladder pressure to calculate detrusor (bladder muscle) pressure. Detrusor pressure is then plotted against volume.

Both the SCMG and CCMG testing modules 1400 and 1700 provide a simple, relatively low cost procedure for recording a cystometrogram (CMG). The SCMG and CCMG testing modules are sterile, disposable assemblies that eliminate the need to disinfect equipment prior to use. This, together with a relatively simple set-up and operational procedure by the physician, greatly reduces the time required to obtain the urodynamic data. The SCMG and CCMG testing modules are more comfortable for the female patient and are more cost effective for the physician. The simplicity of the SCMG and CCMG testing modules, and the control device 102 allows operation with minimal training. Further, when combined in operational use with the SUI testing module 1000, these modules provide a near complete urodynamic diagnostic tool for the physician.

UROFLOMETRY

A uroflometry testing module 2400 can also be removably coupled to control device 102. The module housing of the uroflometry testing module 2400 may be in the form of a plastic disposable cartridge. As shown in Figures 24 and 25, the Uroflometry testing module 2400 includes a single lumen tubing member 2402 having a proximal end 2404 and a distal end 2406 and a channel extending substantially therethrough. A balloon 2408 or other suitable elastomeric element is coupled to the distal end 2406, however, so that the channel of the single lumen tubing member is not open at the distal end. A pressure cushion may also be used in place of the balloon. A collection bucket 2410 is positioned on top of the balloon. The inner surface of the collection bucket may also contain a urinalysis strip which, when wetted by the voided urine, allows for quantitative assessment of standard urinalysis parameters

The diagnostic system including the Uroflometry testing module is operated as follows. The collection bucket is positioned under a commode 2412 to collect urine as the patient voids. Balloon is positioned relative to the bucket so that it substantially supports the bucket. As the bucket fills the pressure in the balloon rises proportionately to the weight of the fluid. When the testing module is coupled to the control device, the proximal end 2404 of the single lumen tubing member 2402 contacts the pressure transducer 128 of the control device 102 so that the pressure within the balloon can be captured and interpreted by the control device. The pressure

data is used to calculate the weight and volume of the fluid (known fluid density). The stored information is then available for review on the display screen, or by a printout through a charging cradle (printer assembly), or downloaded to a PC via a software interface in the charging cradle. Once the test has been completed, the operator
5 disengages the input pendant switches 124, and the urine and collection bucket are discarded.

Operation of the Uroflometry module software subroutine is illustrated in Figures 8a-b. Following module detection 802 and a command to execute the UroFlow Module Subroutine 804, the UroFlow module subroutine begins. The operator
10 is prompted to Enter UroFlow Patient Data 840 necessary for the UroFlow test routine. Once the patient data is collected a UroFlow Scale Zeroing Procedure 841 runs. The operator then enters information necessary to initiate the UroFlow test (UroFlow Test I/O) and the test is started 842. Following the test the software then exits the UroFlow test subroutine and stores the data collected in the Data Storage routine.

15

VAGINAL SPECULUM

Figures 26-28 illustrate a vaginal speculum assembly 2600 for use in the reduction of vaginal prolapse when performing female urodynamic testing, as
20 previously discussed. Uterine or vaginal prolapse occurs when the uterus or pelvic organs drop or become displaced because of weakened pelvic muscles. Prolapse must be reduced to effectively perform urodynamic tests to ensure that no underlying stress urinary incontinence symptoms are masked by the pressure of the vaginal prolapse, which may cause distortion or kinking of the urethral canal. The vaginal speculum
25 assembly 2600 will permit the clinician or physician to perform a urodynamic test procedure with one hand while still reducing vaginal prolapse, as well as properly position the meatus plug device or other catheter within the urethral canal. This prolapse maneuver using the vaginal speculum assembly 2600 during urodynamic testing is especially important prior to surgical repair of the vaginal prolapse, as an
30 undiagnosed case of stress urinary incontinence may surface following prolapse surgery. The urodynamic testing being performed using the vaginal speculum assembly in this manner allows the surgeon to determine if additional stress urinary incontinence (SUI) surgery should be performed at the time of prolapse repair.

Current medical practice calls for the use of a vaginal speculum secured in place in order to reduce the prolapse. For example, U.S. Patent Nos. 5,997,474 and 6,048,308 describe specula specifically designed for vaginal examination and treatment. U.S. Patent No. 6,120,438 discloses a vaginal retractor device designed to hold back the vaginal wall during an exam or surgical procedure. Often, surgical tape is necessary to hold the speculum in place, as the physician's hands cannot hold the speculum in place while performing a particular urodynamic procedure. None of the prior art speculum devices integrate the use of urodynamic equipment.

With reference to Figures 26 and 27, the vaginal speculum assembly 2600 includes a connector member 2602 for coupling an insertion device assembly, such as a meatus plug device 1102, or catheter 1422 and related elements to the vaginal speculum. The vaginal speculum can be of any type well known in the art. In the illustrated embodiment, the vaginal speculum includes an upper arm 2604, a lower arm 2606, and a hinge member 2608 for joining the upper and lower arms together. The vaginal speculum also includes a handled member 2610 being integrally attached, and preferably substantially perpendicular aligned to the lower arm. The vaginal speculum 2600 further includes a locking bar device 2612 connected to the upper arm 2606 for locking the upper and lower arms in an open position, as shown in Figure 28. The upper arm 2604 includes a posterior end 2614 with a pair of arm mounting openings 2616 therein. The connector member 2602 includes a flexible band 2618. The flexible band at one end 2620 includes a pair of mounting openings 2622 and at the other end 2624 a connector element 2626. The mounting openings 2622 of the flexible band 2618 are aligned with the arm mounting openings 2616 of the upper arm for receiving a pair of mounting screws 2628 therein in order to attach the connector member 2602 to the vaginal speculum 2600. During use, the connector element can be coupled to the meatus plug device or catheter as shown in Figure 26.

Although a particular embodiment of the connector member 2602 is illustrated and described herein, those skilled in the art will recognize that various other embodiments are also possible to provide a means by which to removably couple a device that is inserted into the urethral canal to the speculum so as to hold it in place within the patient.

In operation, the vaginal speculum assembly 2600 can be cooperatively used in conjunction with the urodynamic system disclosed herein. For example, it may

be used in conjunction with a urodynamic system including a SUI testing module 1000 in the performance of the urodynamic testing procedure for stress urinary incontinence (SUI), such as the measuring of urethral resistance pressure (URP) as previously described. In reference to Figure 28, the physician positions the vaginal speculum assembly 2600, such that it is fully inserted within vaginal canal 2650 wherein the upper and lower arms 2604, 2606 are fully opened and pressed against the vaginal walls 2650w for reducing the patient's vaginal prolapse. The physician then locks the upper and lower arms of the vaginal speculum in the fully opened configuration (see Figure 28) via the locking bar device 2612, and adjusts the connector member 2602 so that the insert member will be aligned with the urethral canal. The remaining operational steps are exactly the same as the operational steps described above in connection with individual testing modules.

Although the portable medical system disclosed herein has been described in conjunction with diagnostic testing, it is to be understood that the system can also be used in conjunction with therapies and/or surgical procedures for treating urinary incontinence, such as placement of a sling, placement of bulking agents, shrinkage of tissue etc. In this regard, the testing described herein can be used before, during and/or after these procedures to ensure success of the procedures, for example, to ensure correct placement and/or tensioning of a sling.

Although exemplary embodiments and methods for use have been described in detail above, those skilled in the art will understand that many variations are possible without departing from the spirit and scope of the invention, which is limited only by the appended claims.

What is claimed is:

1. A vaginal speculum assembly comprising:

5 an upper arm and a lower arm movably coupled to one another by a coupling member positioned at a posterior end, each of the upper and lower arms being dimensioned for at least partial insertion into a patient's vagina, and when so inserted, being dimensioned to exert pressure on at least upper and lower portions of an inner surface of the patient's vaginal wall, the coupling element permitting movement of the upper and lower arms relative to one another to adjust their positioning relative to the vaginal wall;

10 a locking device coupled to at least one of the upper and lower arms that, when activated, substantially prevents movement of the upper and lower arms relative to one another; and

15 a connector member extending outwardly from a posterior end of the upper arm, the connector member being configured to engage a predetermined portion of an insert device assembly including an insert member to be inserted into the patient's urethral canal to thereby maintain the insert device assembly in position relative to the vaginal speculum assembly.

20 2. The vaginal speculum assembly according to claim 1, wherein at least a portion of the connector member is comprised of a flexible material.

25 3. The vaginal speculum assembly according to claim 2, wherein the flexible portion can be manually manipulated to engage and grasp the insert device assembly.

30 4. The vaginal speculum assembly according to claim 1, wherein the coupling member is a hinge member, and said relative movement is pivotal movement about said hinge member.

5. The vaginal speculum assembly according to claim 1, wherein the coupling element permits movement of at least a portion of the upper and lower arms relative to one another in a substantially radial direction.
- 5 6. The vaginal speculum assembly according to claim 1, wherein, when the upper and lower arms are inserted into the patient's vagina, the connector member is further configured to maintain the insert device assembly in position relative to the vaginal speculum assembly so that the insert member extends into the patient's urethral canal.
- 10 7. The vaginal speculum assembly according to claim 6, wherein the insert member is a meatus plug member.
8. The vaginal speculum assembly according to claim 6, wherein the insert member is a catheter.
- 15 9. A method for assessing urinary function comprising the steps of:
- positioning a vaginal speculum assembly relative to a patient's vagina so that upper and lower arms of the vaginal speculum assembly support at least upper and lower portions of an inner wall of the patient's vagina;
- 20 locking the upper and lower arms in position relative to one another;
- inserting an insert member of an insert device assembly at least partially into the patient's urethral canal;
- 25 coupling the insert device assembly to the vaginal speculum assembly so that the insert member remains positioned within the patient's urethral canal; and
- 30 performing a urodynamic assessment utilizing the insert device assembly.
10. The method according to claim 9, wherein the urodynamic assessment is a stress urinary incontinence test.

11. The method according to claim 9, wherein the urodynamic assessment is a cystometrogram test.
- 5 12. The method according to claim 9, wherein the insert member is a meatus plug device.
13. The method according to claim 9, wherein the insert member is a catheter.
- 10 14. The method according to claim 9, wherein the step of performing a urodynamic assessment further comprises the steps of infusing fluid into the patient's urinary tract, and measuring pressure within the patient's urinary tract as fluid is infused therein.
- 15 15. The method according to claim 9, wherein the coupling step further comprises the step of coupling a connector member that extends outwardly from a posterior end of the upper arm to the insert device assembly.
16. The method according to claim 15, wherein in the positioning step, the upper and lower arms of the vaginal speculum assembly reduce the patient's vaginal prolapse.
- 20 17. A vaginal speculum assembly comprising:
- an upper arm and a lower arm dimensioned for at least partial insertion into a patient's vagina and movably coupled to one another at a posterior end, the upper and lower arms being movably relative to one another between a first collapsed position and a second open position, wherein when inserted into the patient's vagina and in the open position, the upper and lower arms exert pressure on upper and lower portions of an inner surface of the patient's vaginal wall;
- 25
- a locking device coupled to at least one of the upper and lower arms that, when activated, substantially prevents movement of the upper and lower arms relative to one another; and
- 30

a connector member extending outwardly from a posterior end of the upper arm, the connector member being configured to engage a predetermined portion of an insert device assembly including an insert member to be inserted into the patient's urethral canal to thereby maintain the insert device assembly in position relative to the vaginal speculum assembly.

18. The vaginal speculum assembly according to claim 17, wherein the connector member is further configured to, when the upper and lower arms are inserted into the patient's vagina, to maintain the insert member in position within the patient's urethral canal.


19. The vaginal speculum assembly according to claim 18, wherein the insert member is a meatus plug device.

20. The vaginal speculum assembly according to claim 18, wherein the insert member is a catheter.


21. The vaginal speculum assembly according to claim 17, wherein at least the portion of the connector member configured to engage the insert device assembly is comprised of a flexible material.

22. The vaginal speculum assembly according to claim 17, wherein the coupling member is a hinge member.


1/41

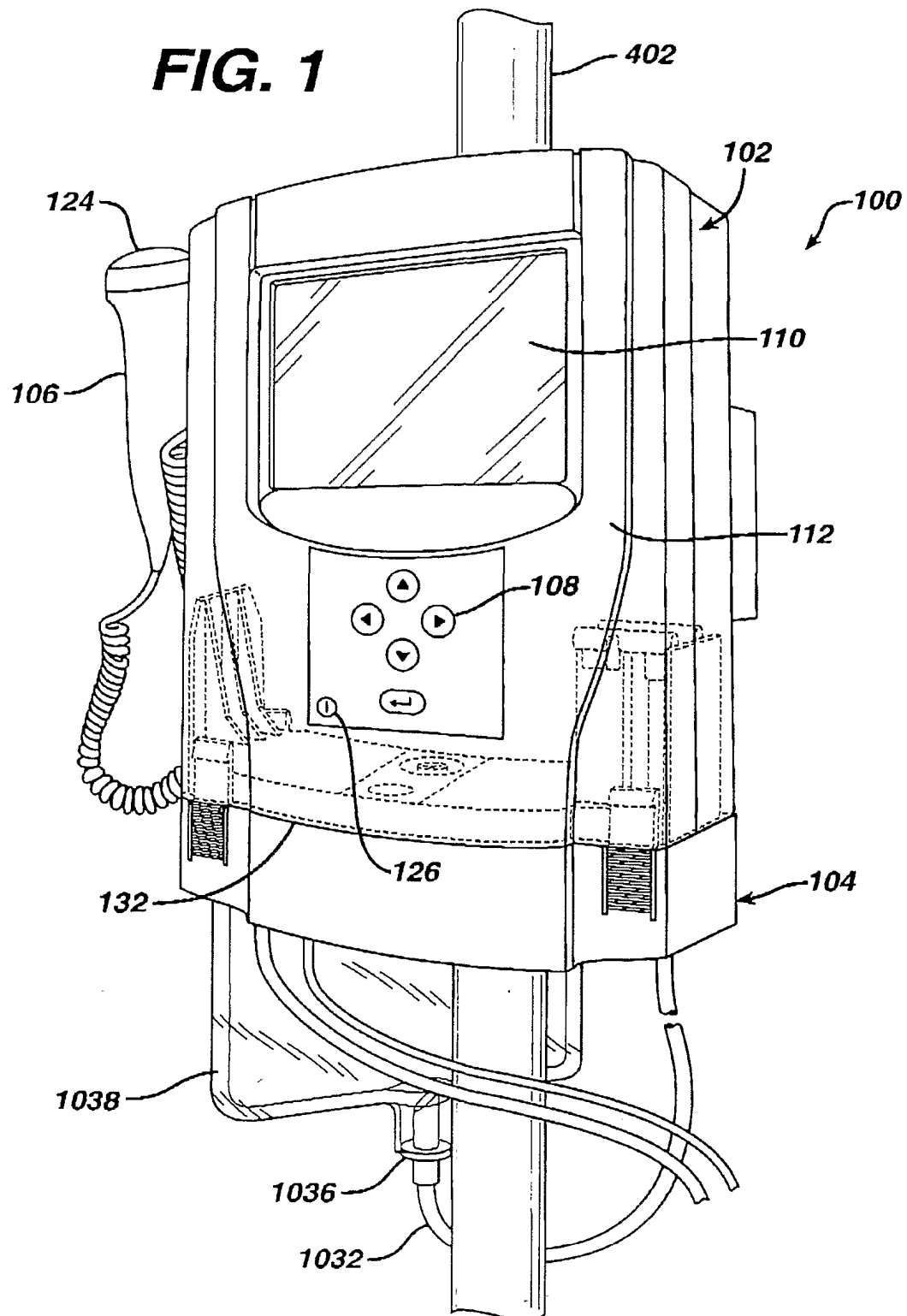FIG. 1


2/41

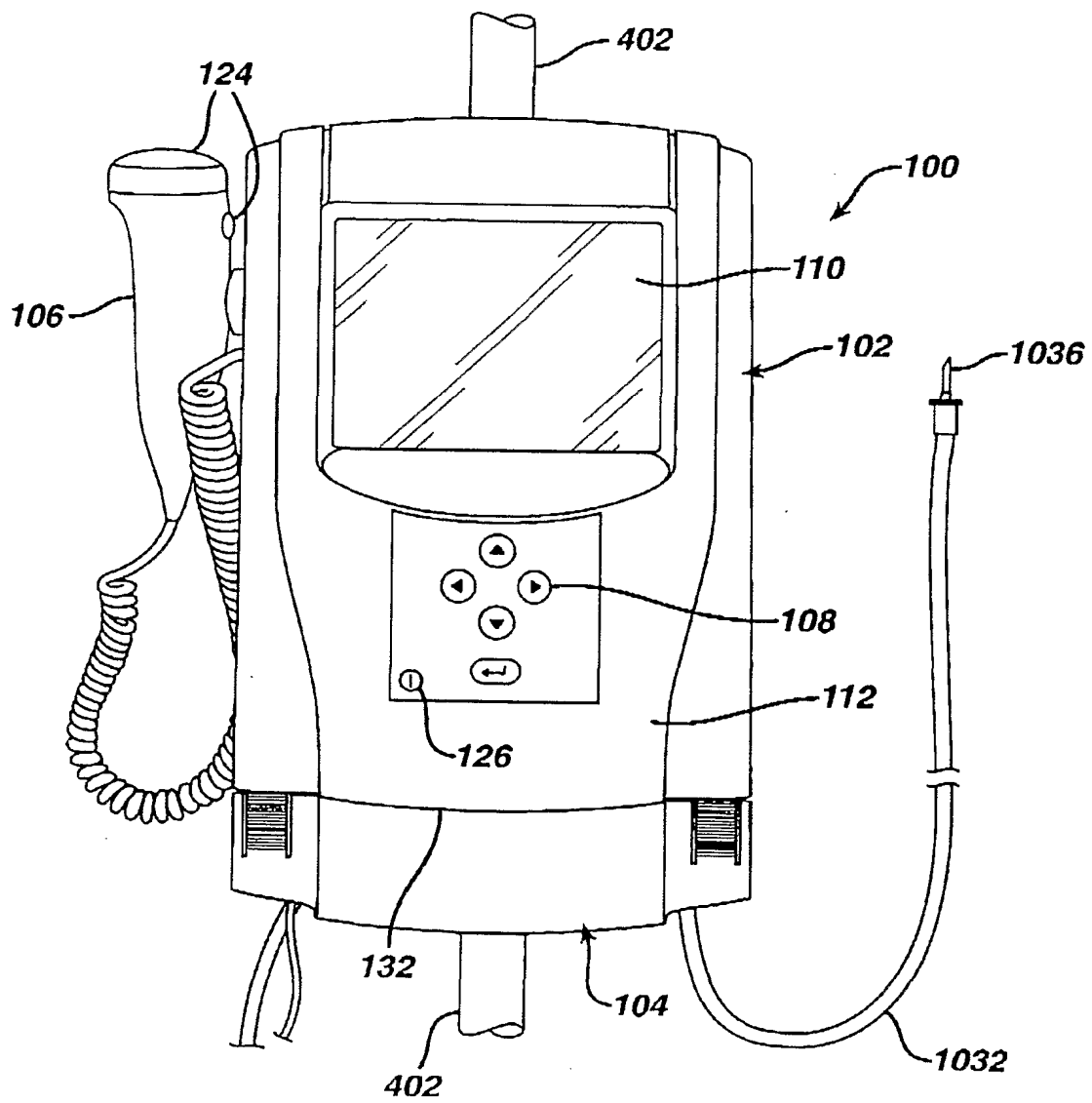FIG. 1a


3/41

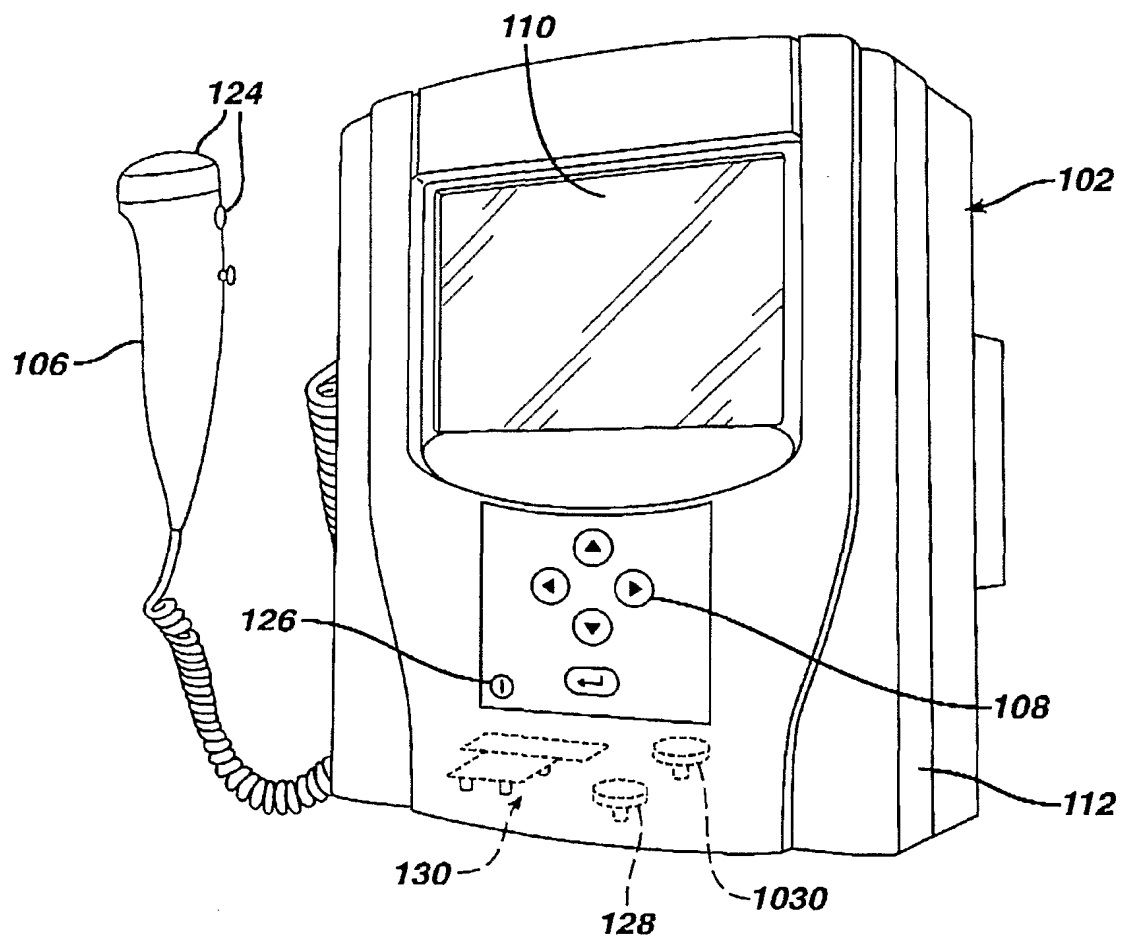FIG. 2


4/41

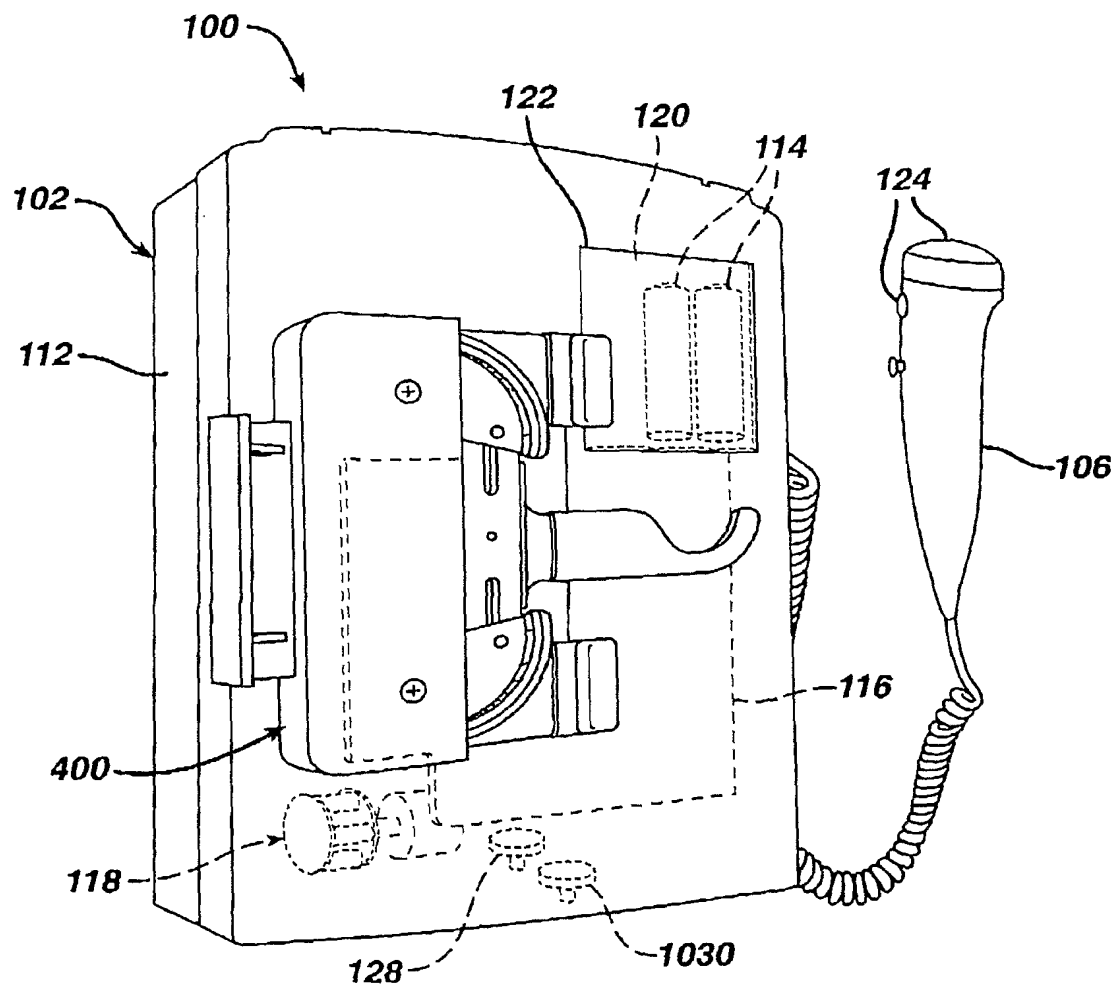FIG. 3


5/41

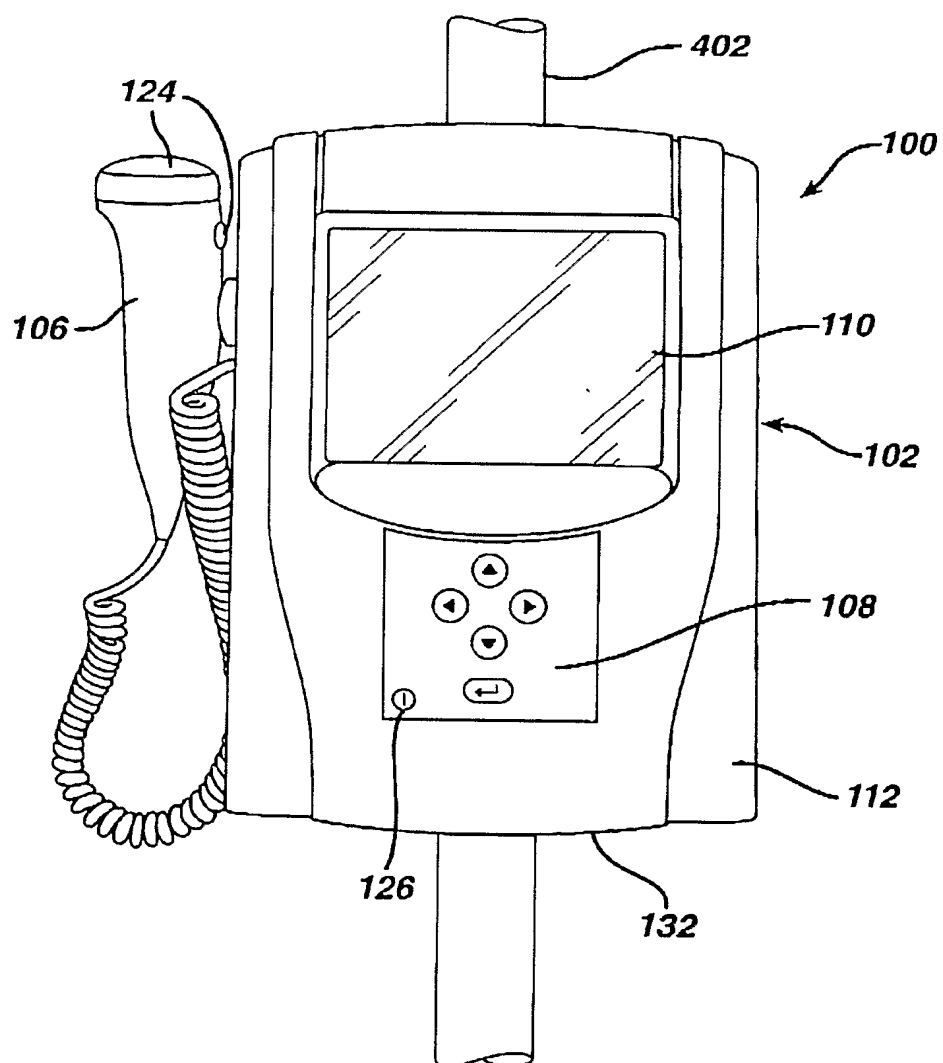FIG. 4

6/41

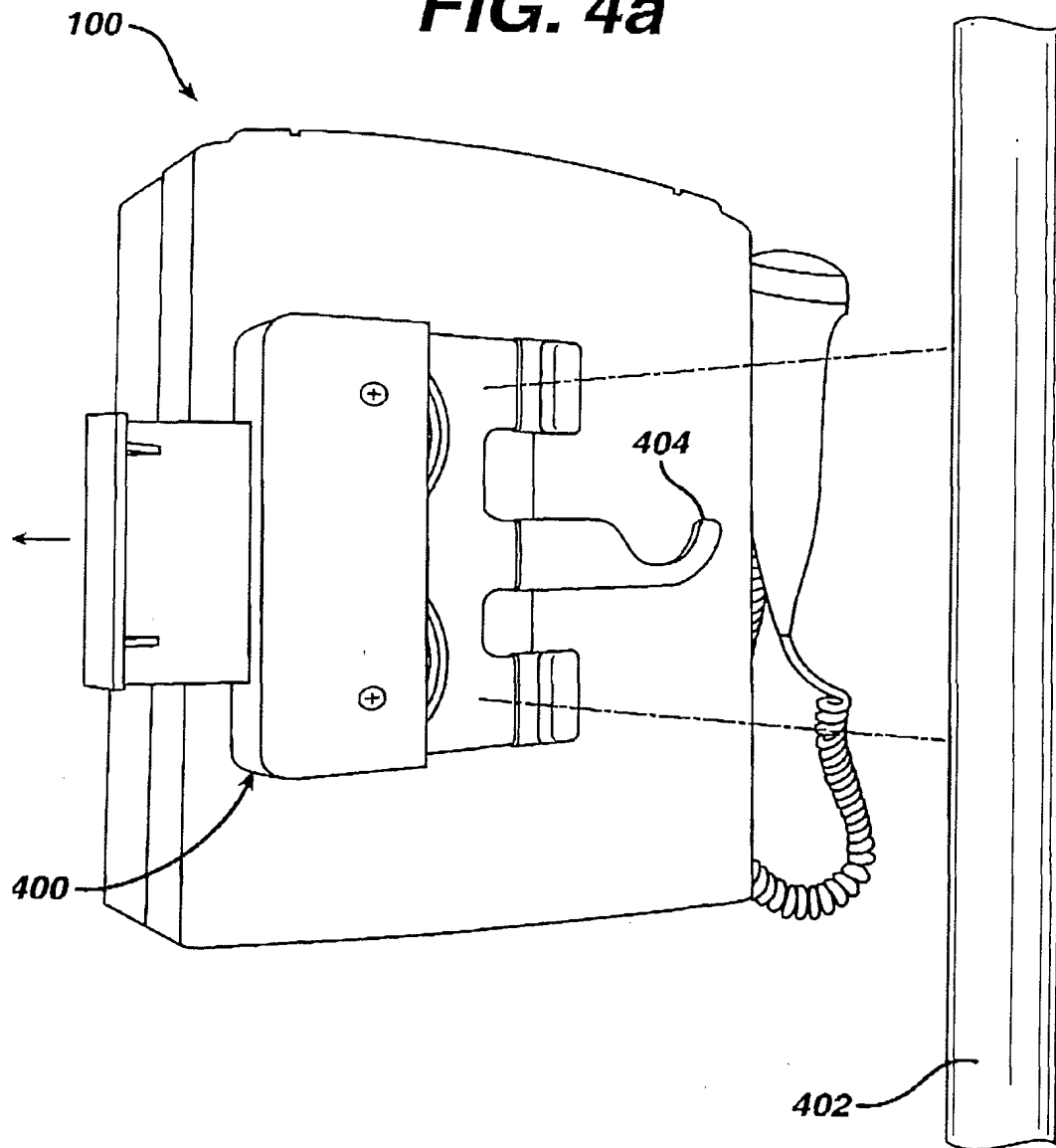FIG. 4a

7/41

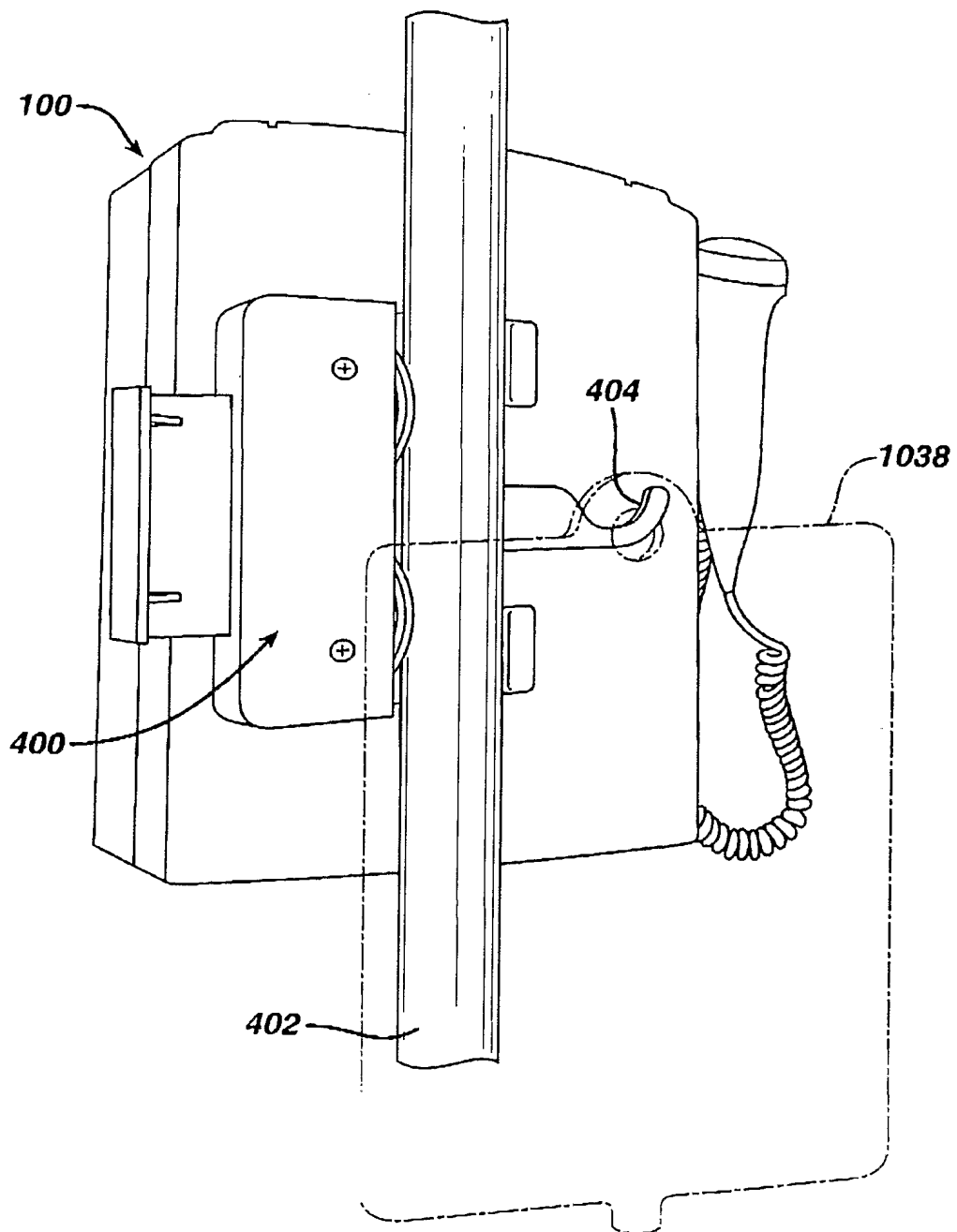FIG. 4b

8/41

FIG. 4c



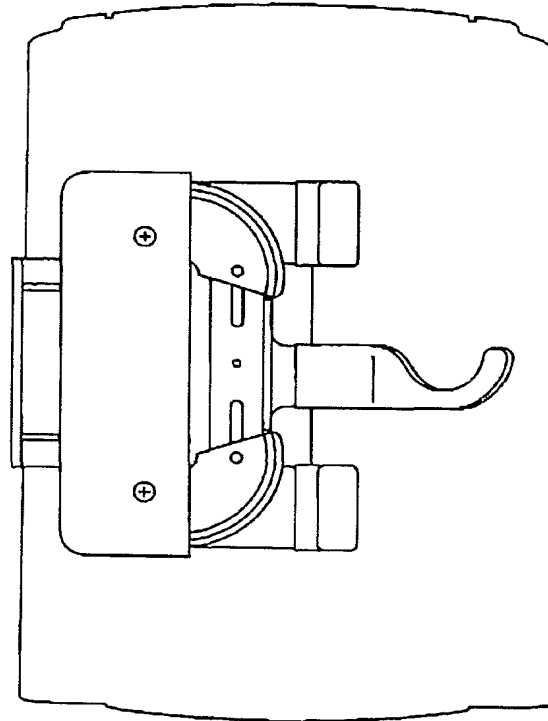


FIG. 4d

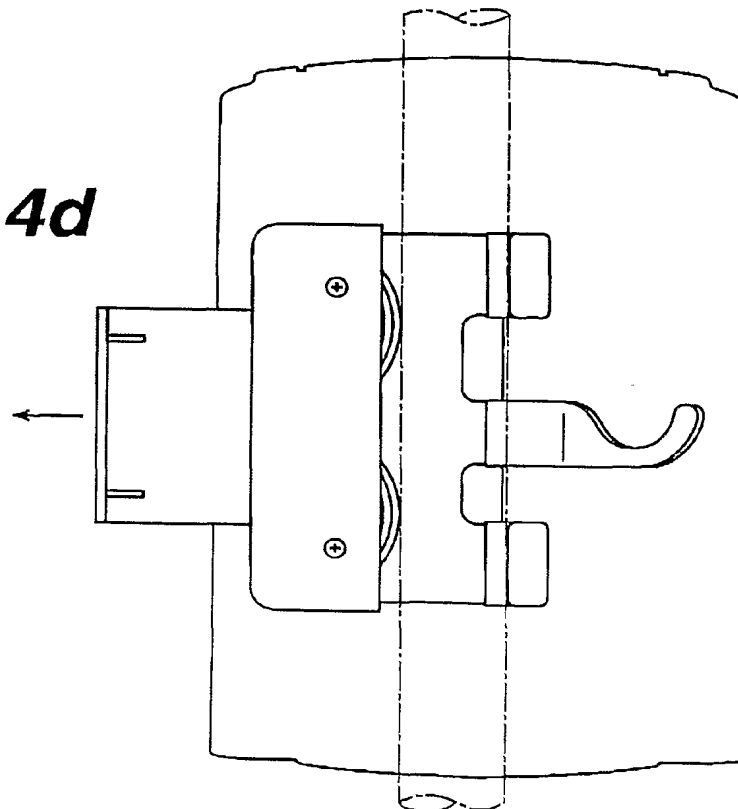

9/41

FIG. 5

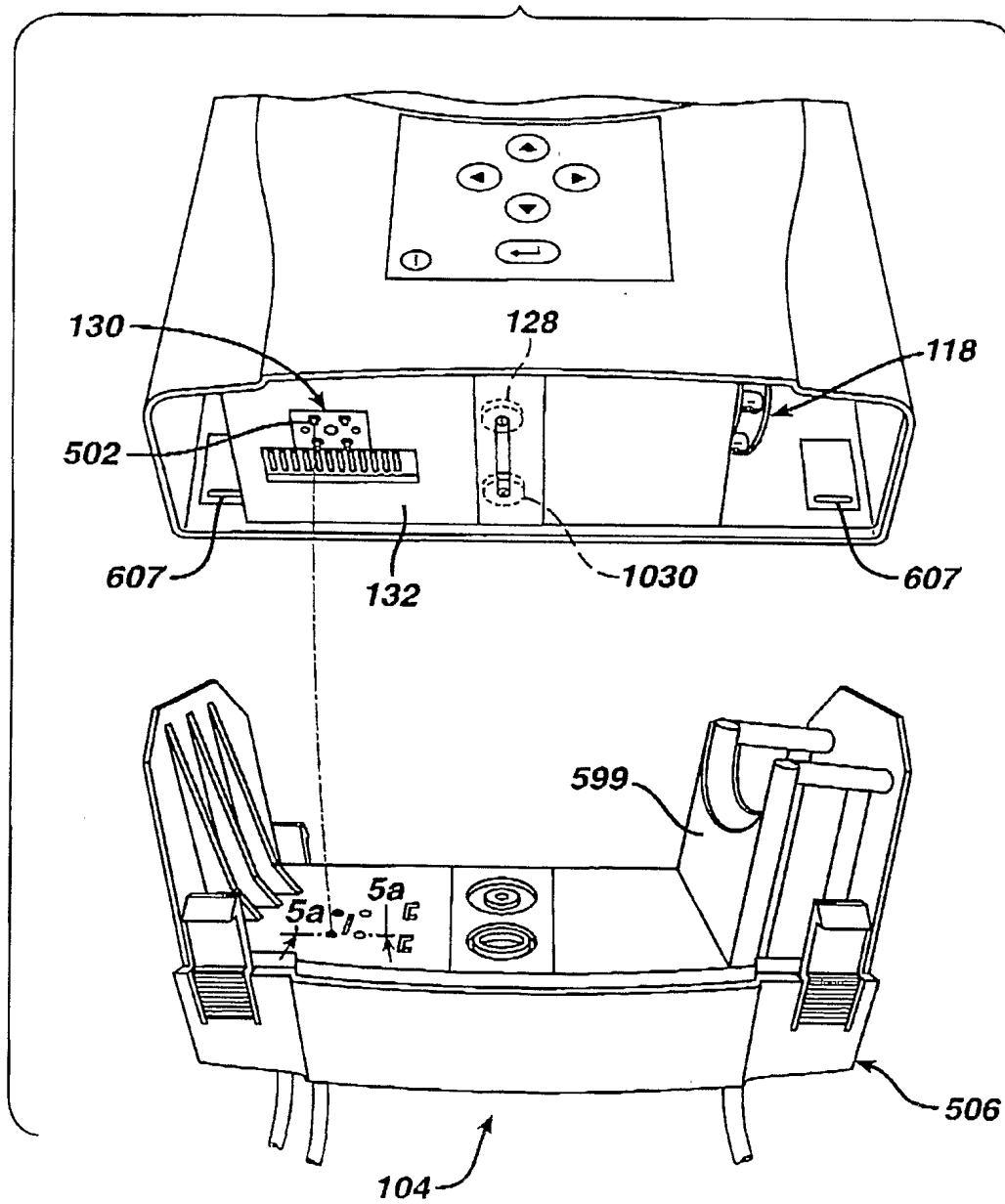

10/41

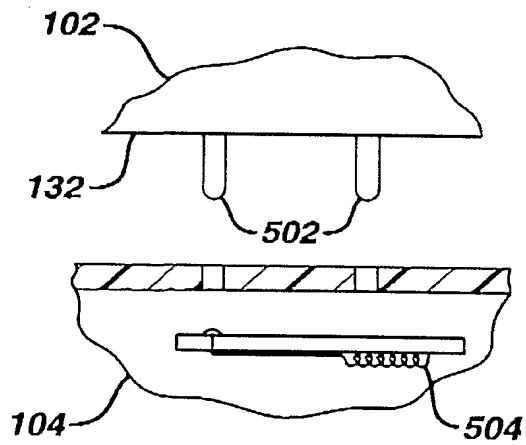


FIG. 5a


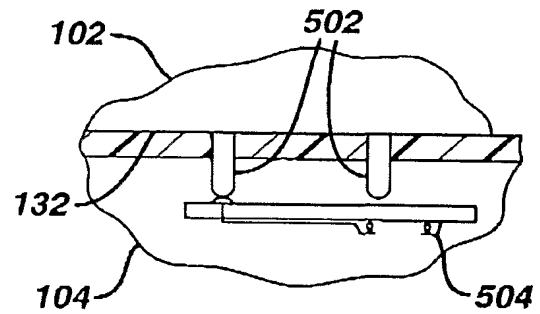


FIG. 5b

FIG. 6

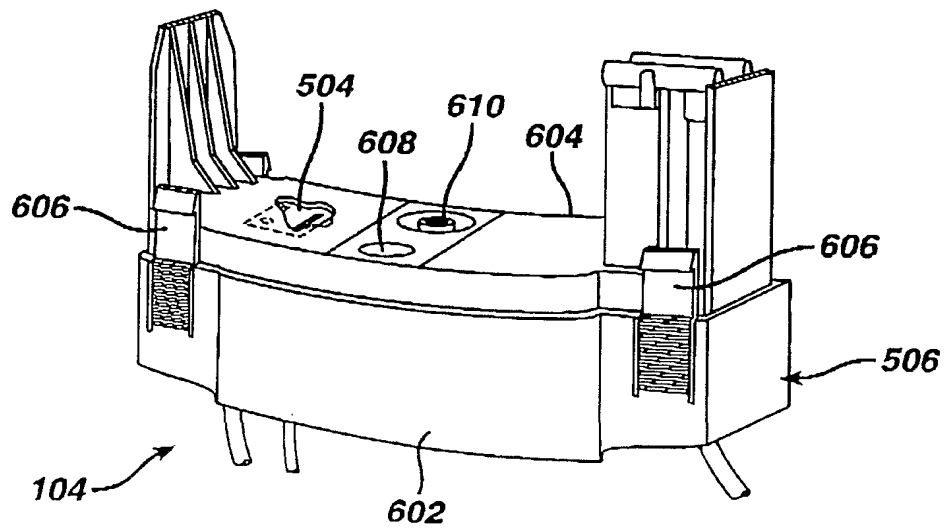

11/41

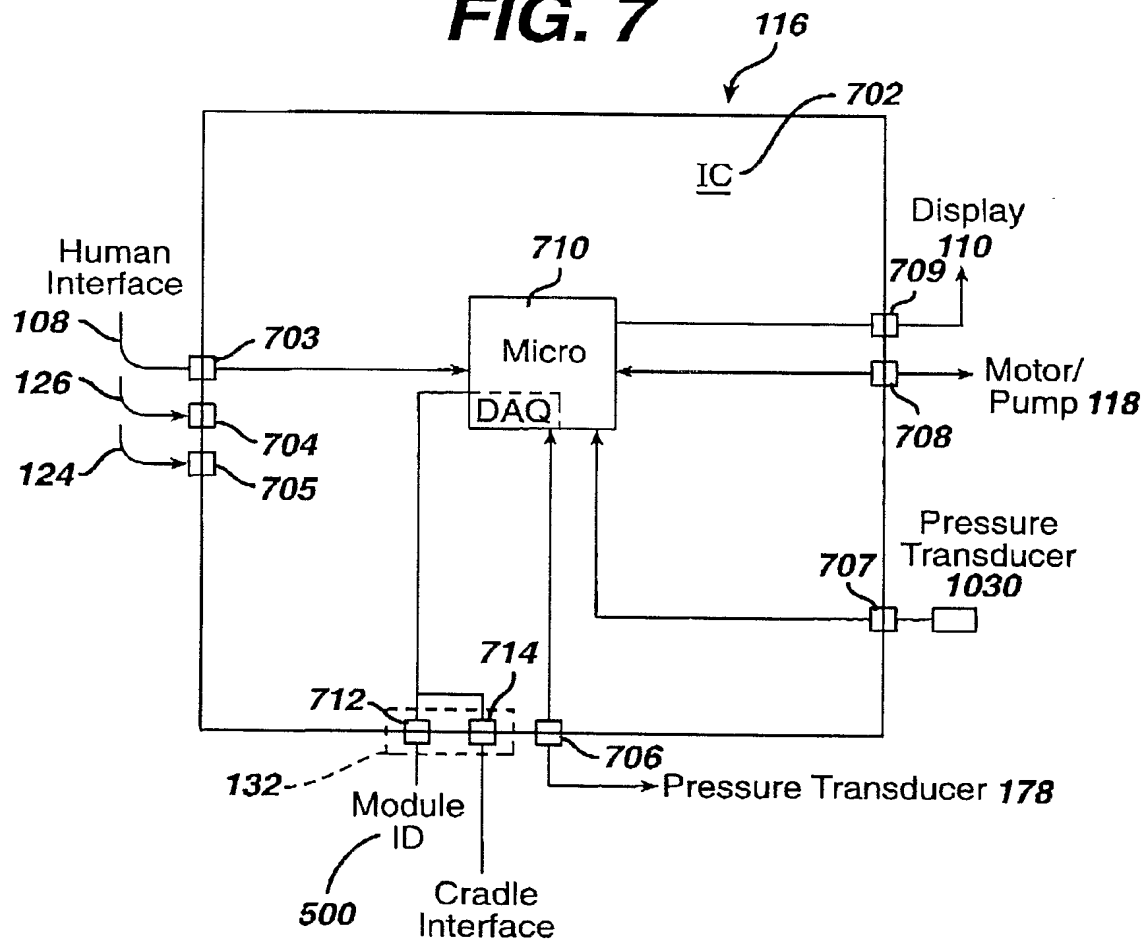
FIG. 7

FIG. 8a


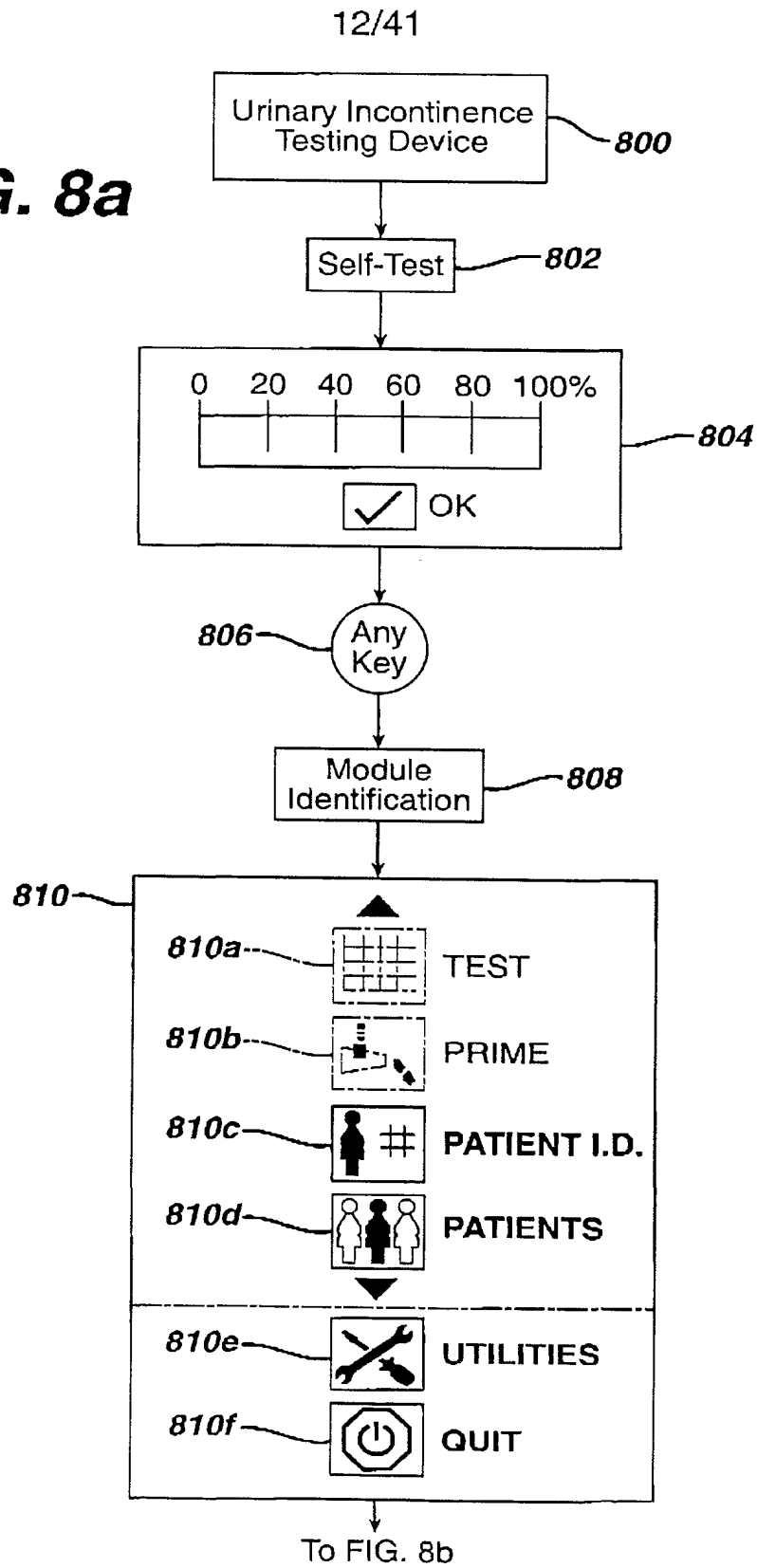


FIG. 8b

14/41




FIG. 8c

FIG. 8d

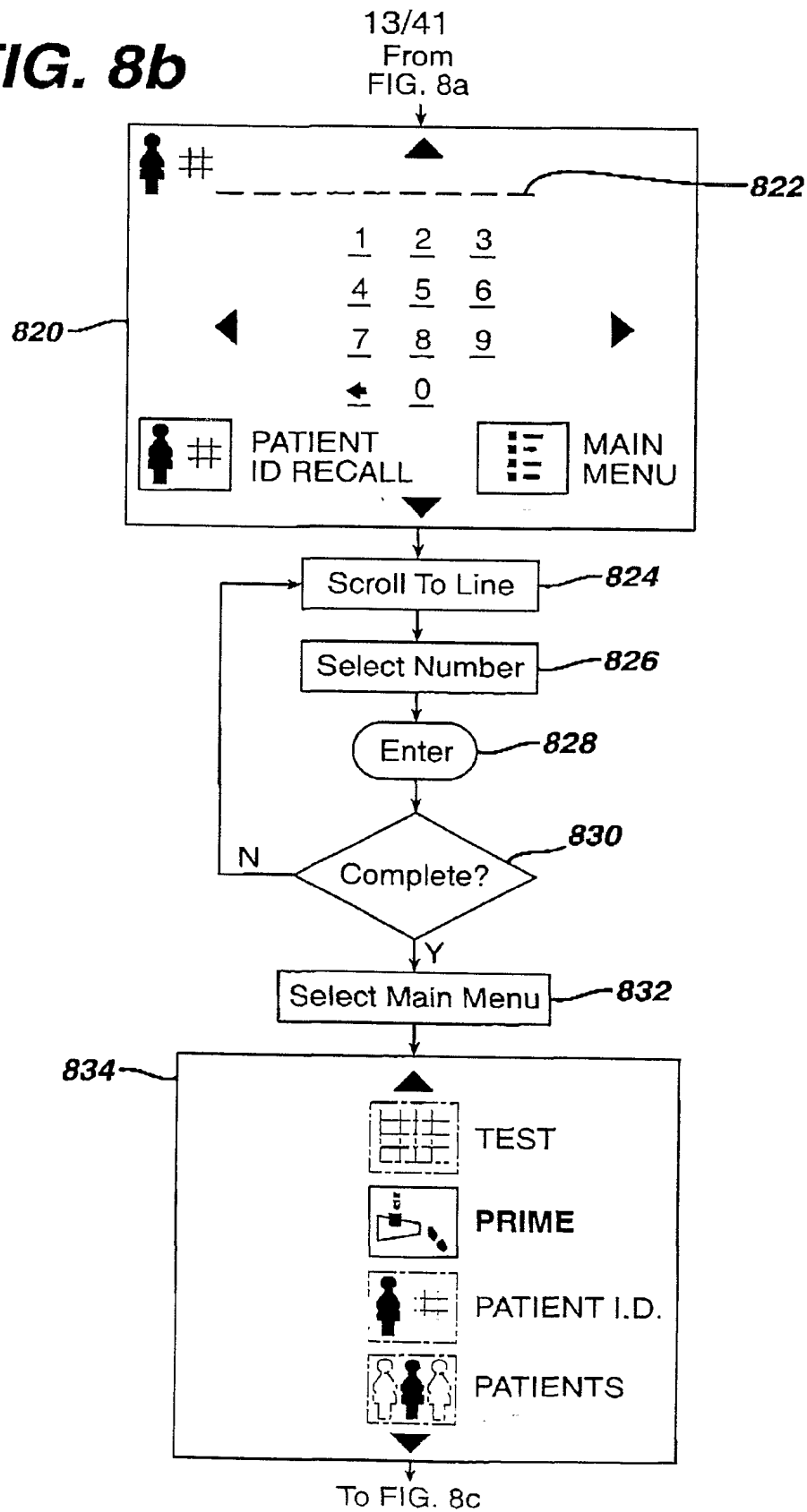

16/41

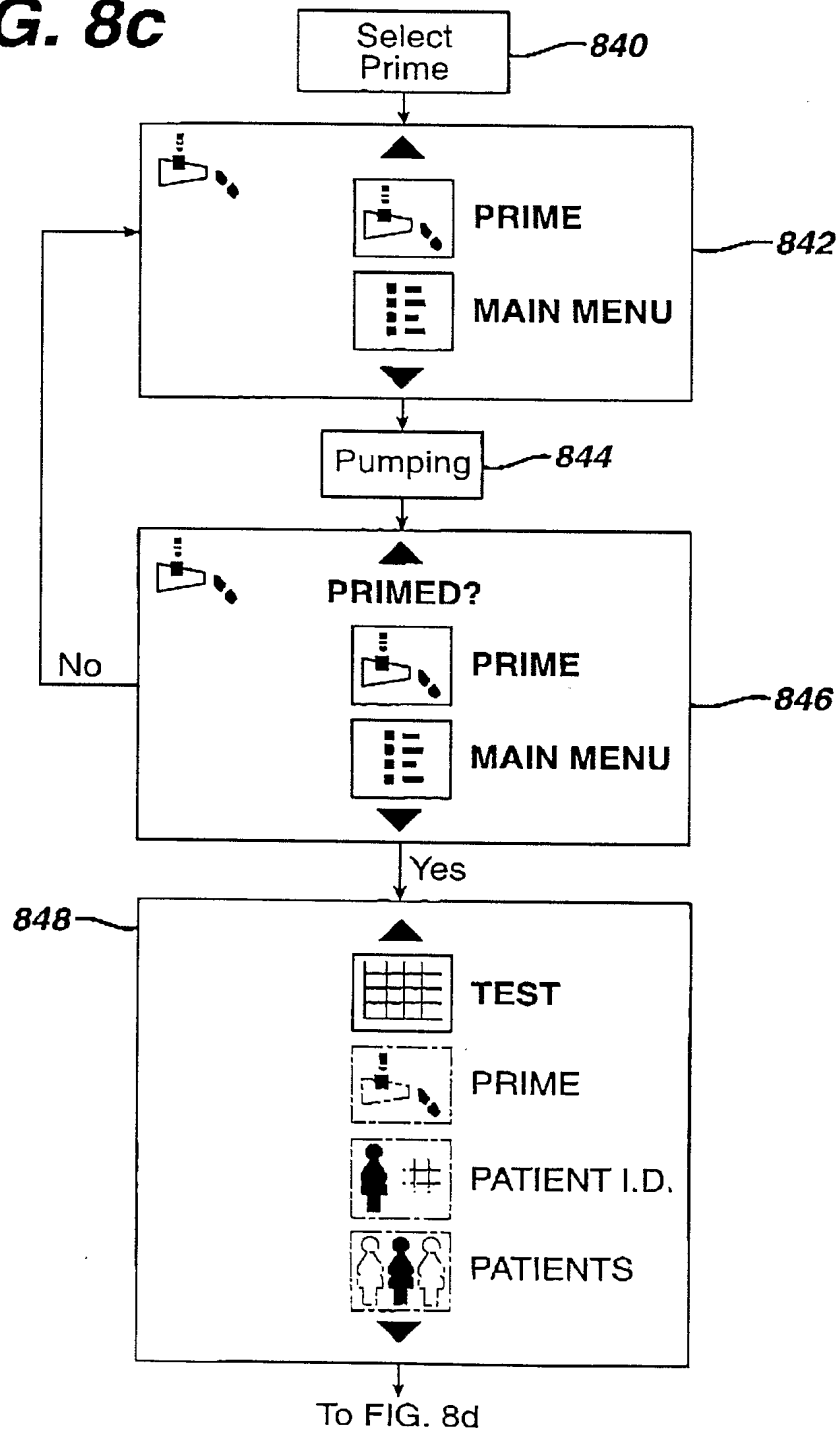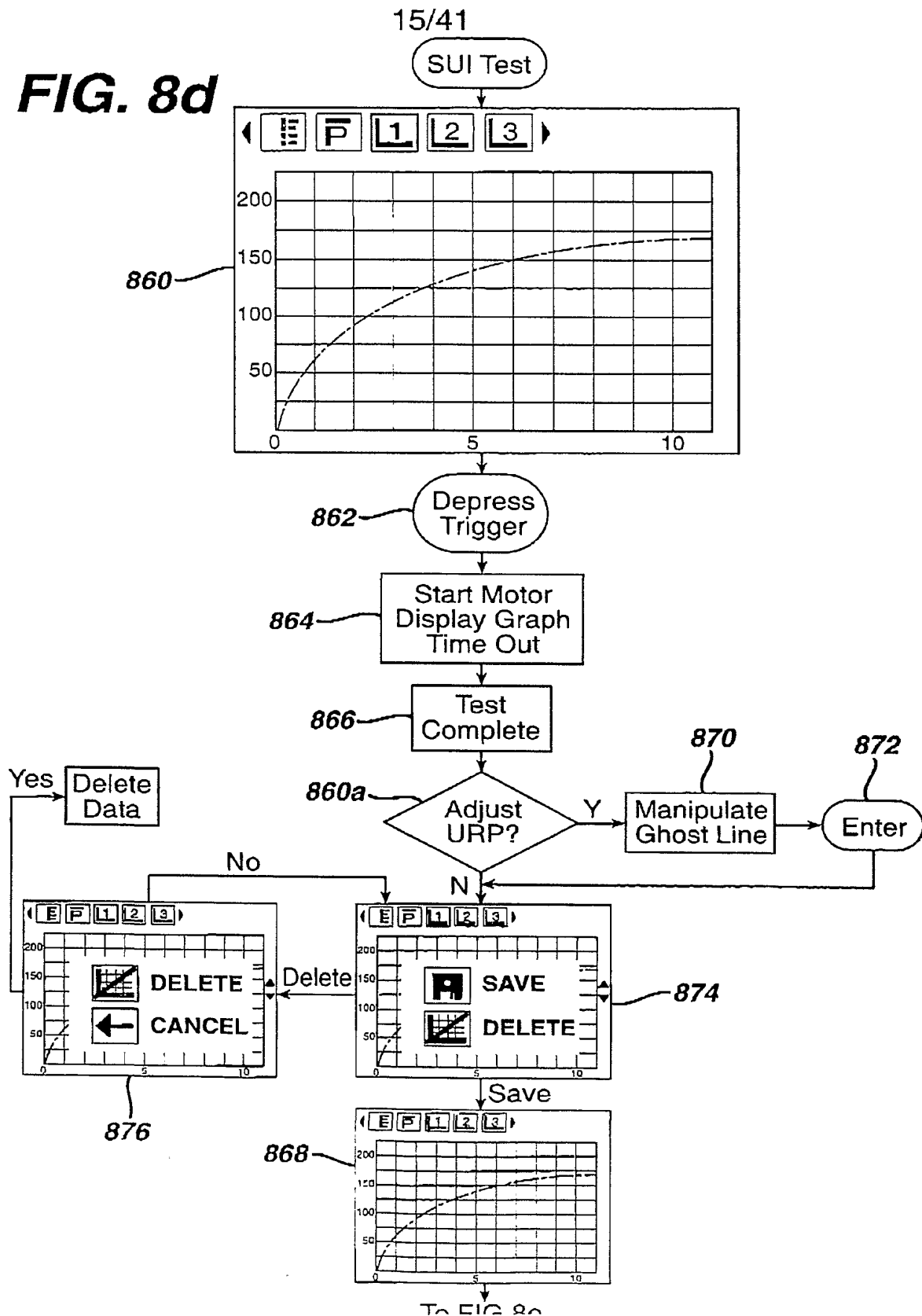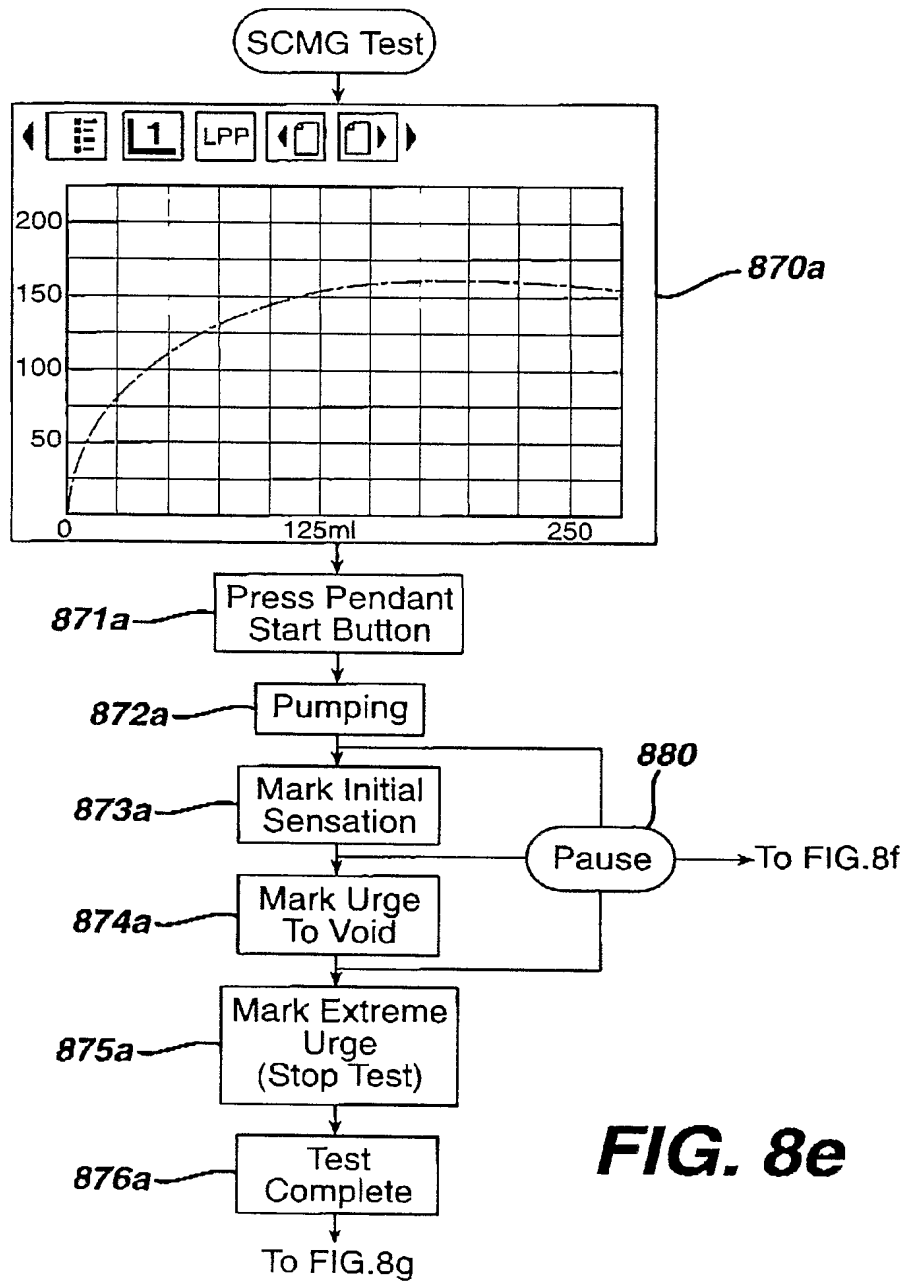

17/41


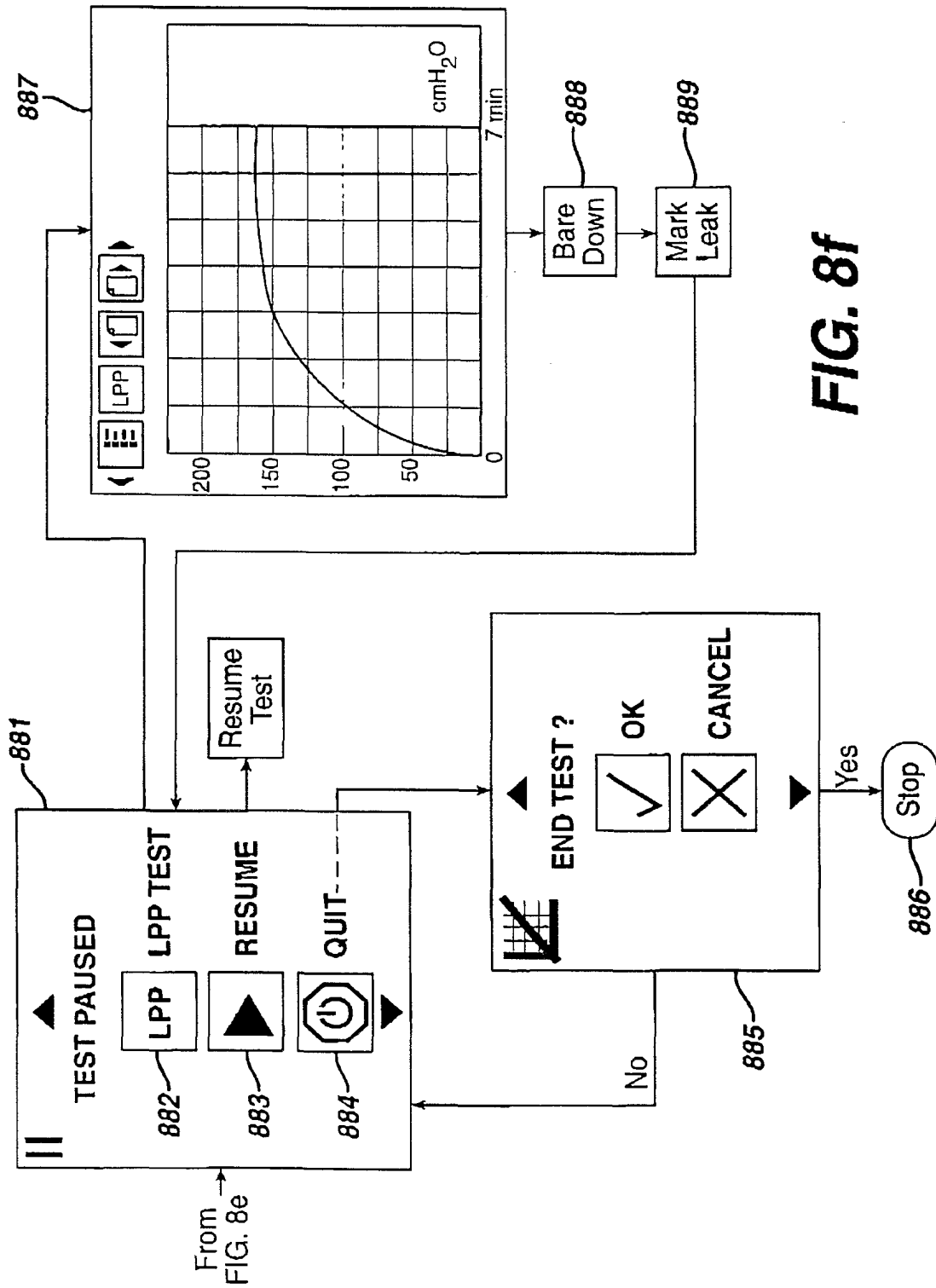


FIG. 8f


18/41

FIG. 8g

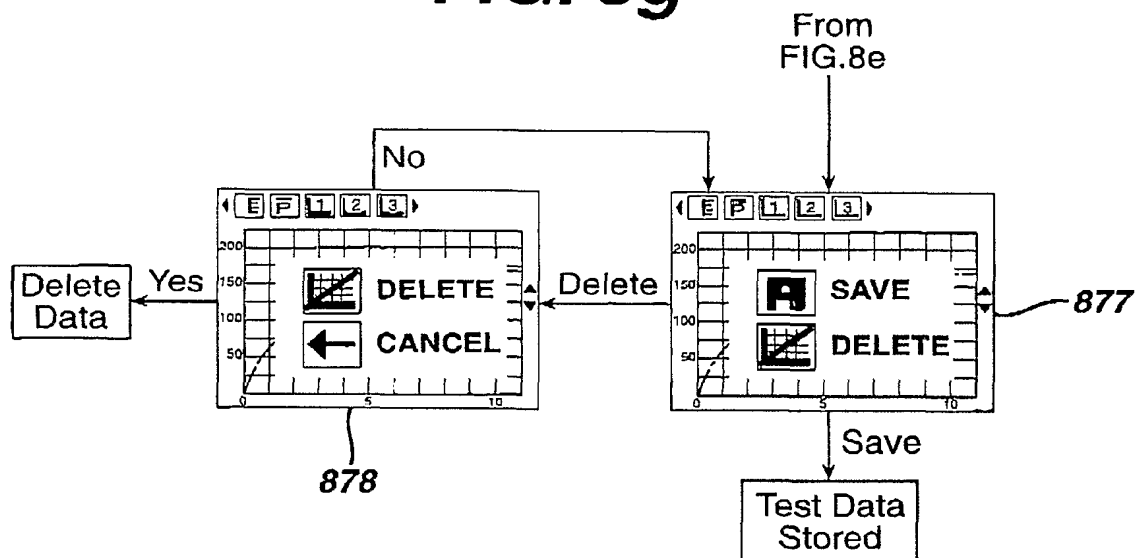


19/41

FIG. 8h

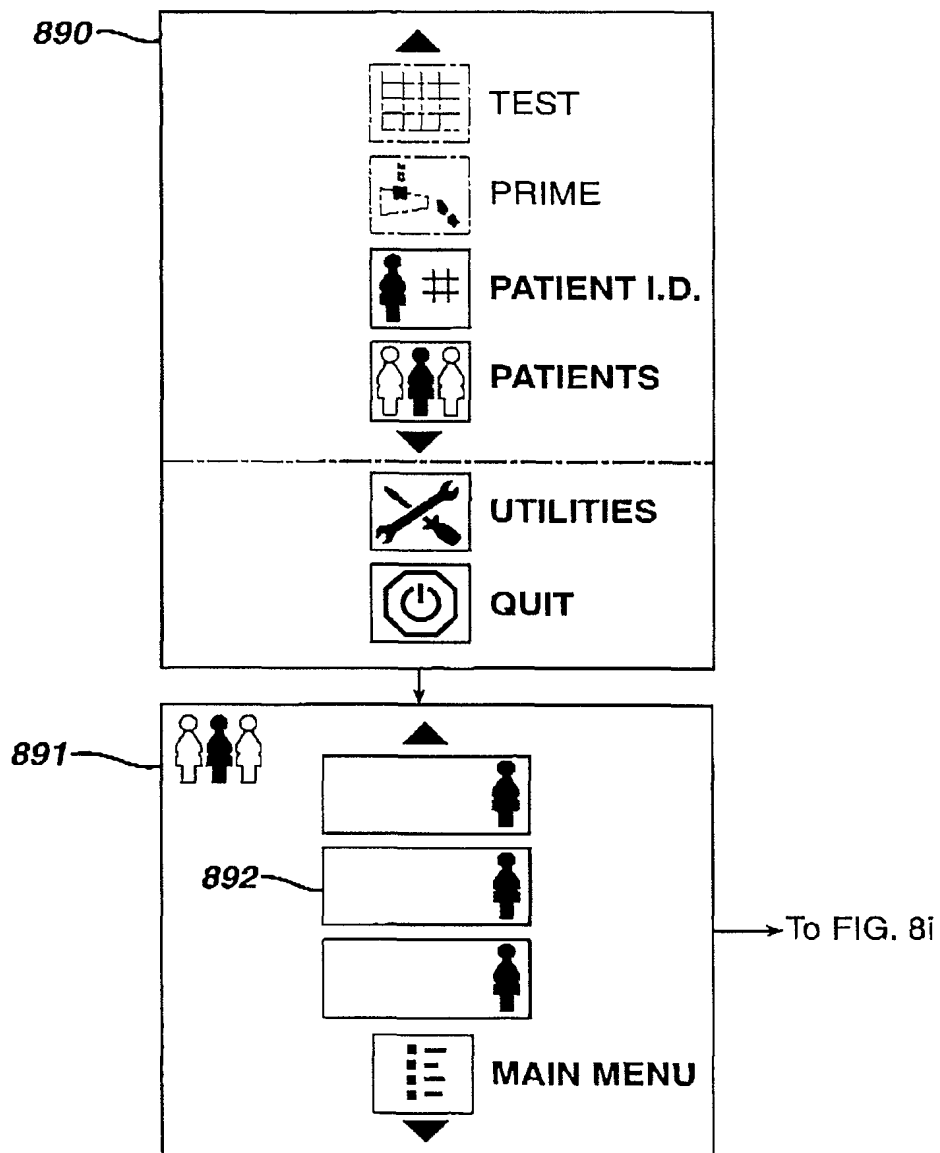


20/41


FIG. 8i

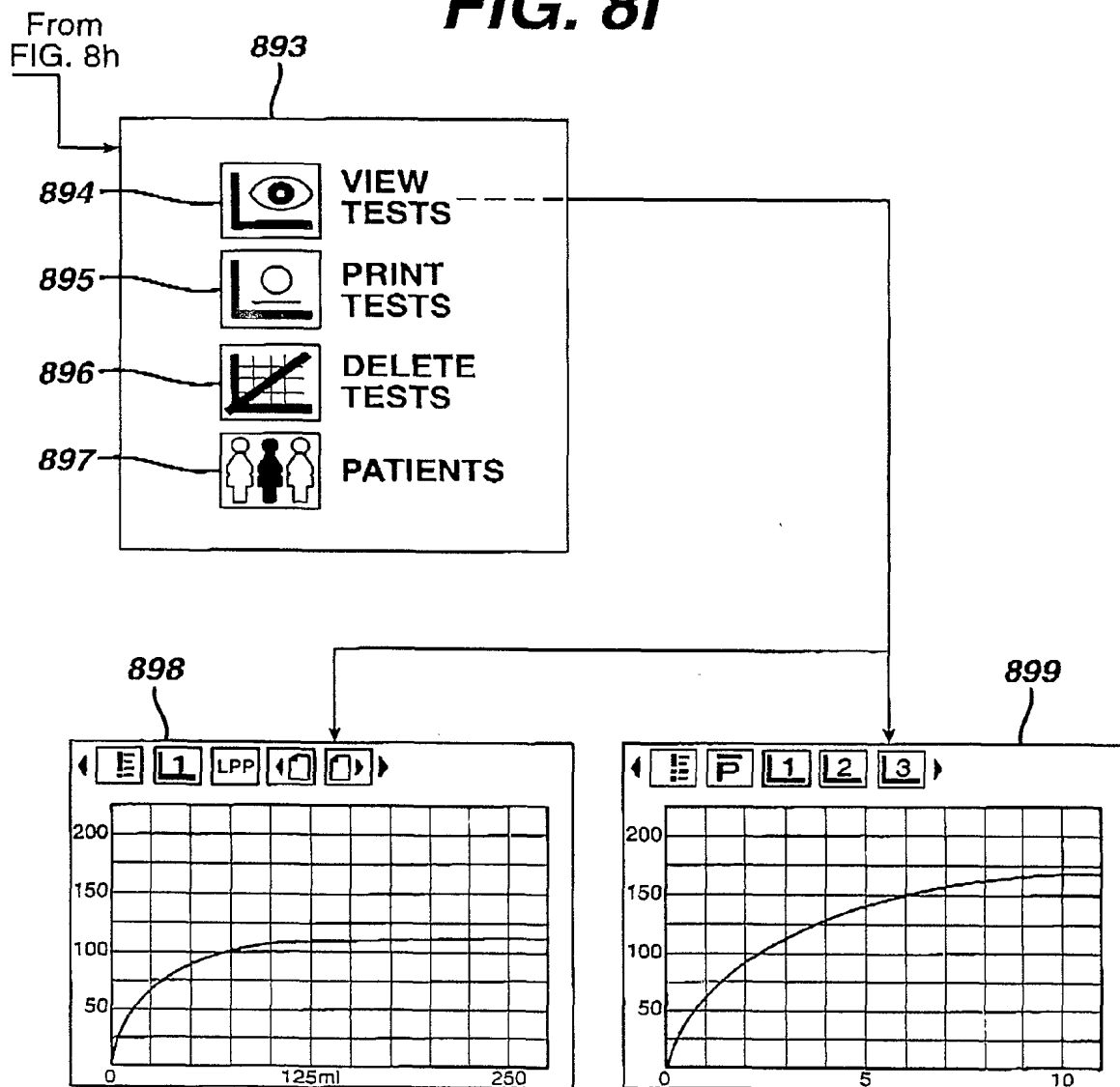

21/41

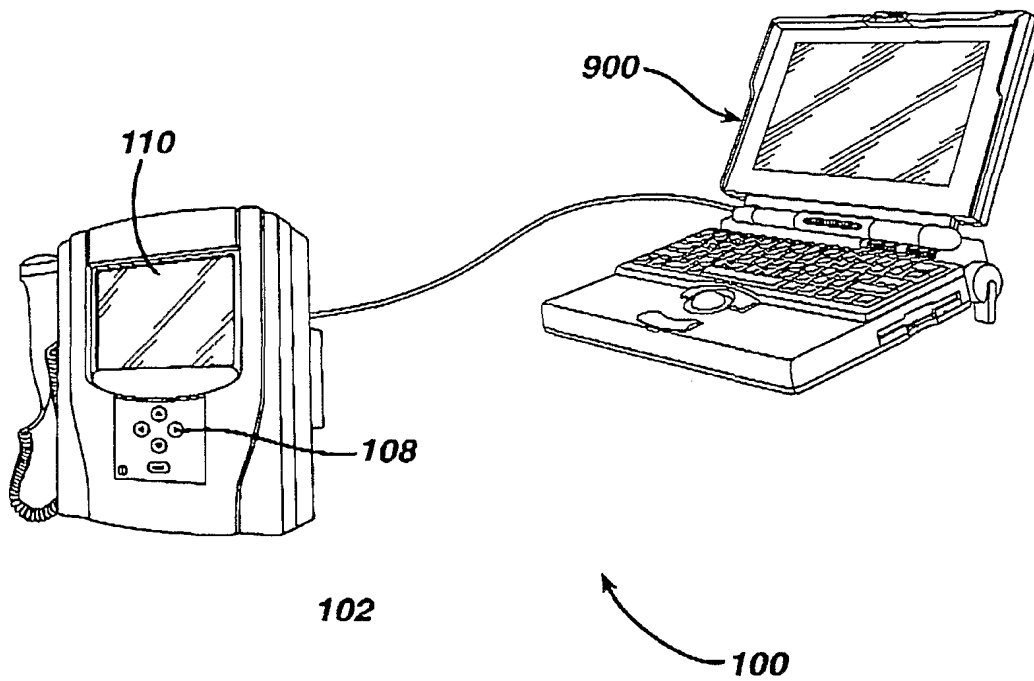FIG. 9

22/41

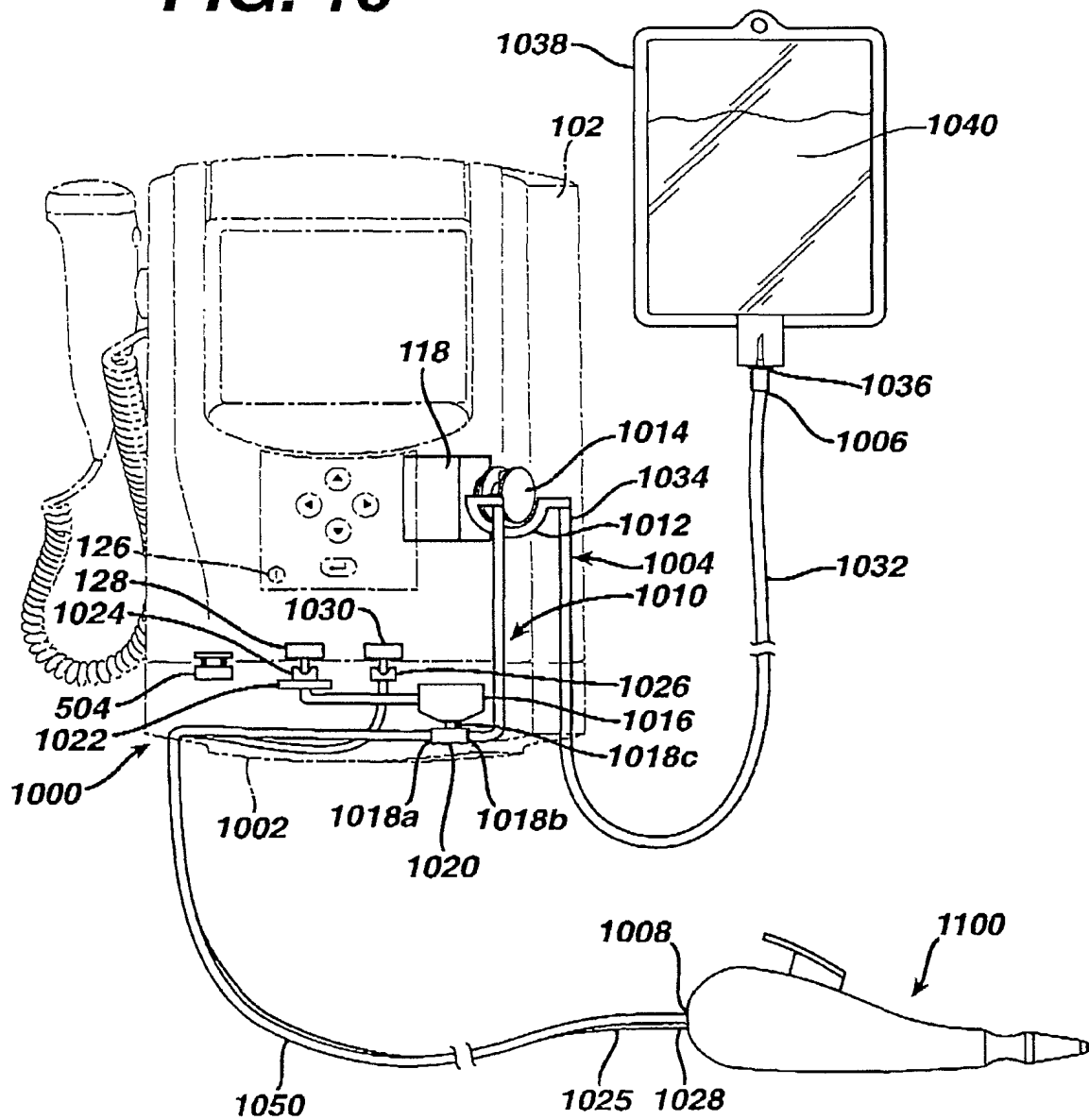FIG. 10

23/41

FIG. 10a







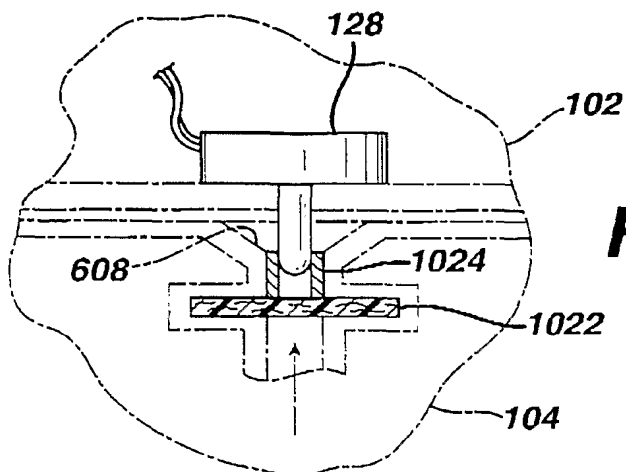


FIG. 10b

24/41

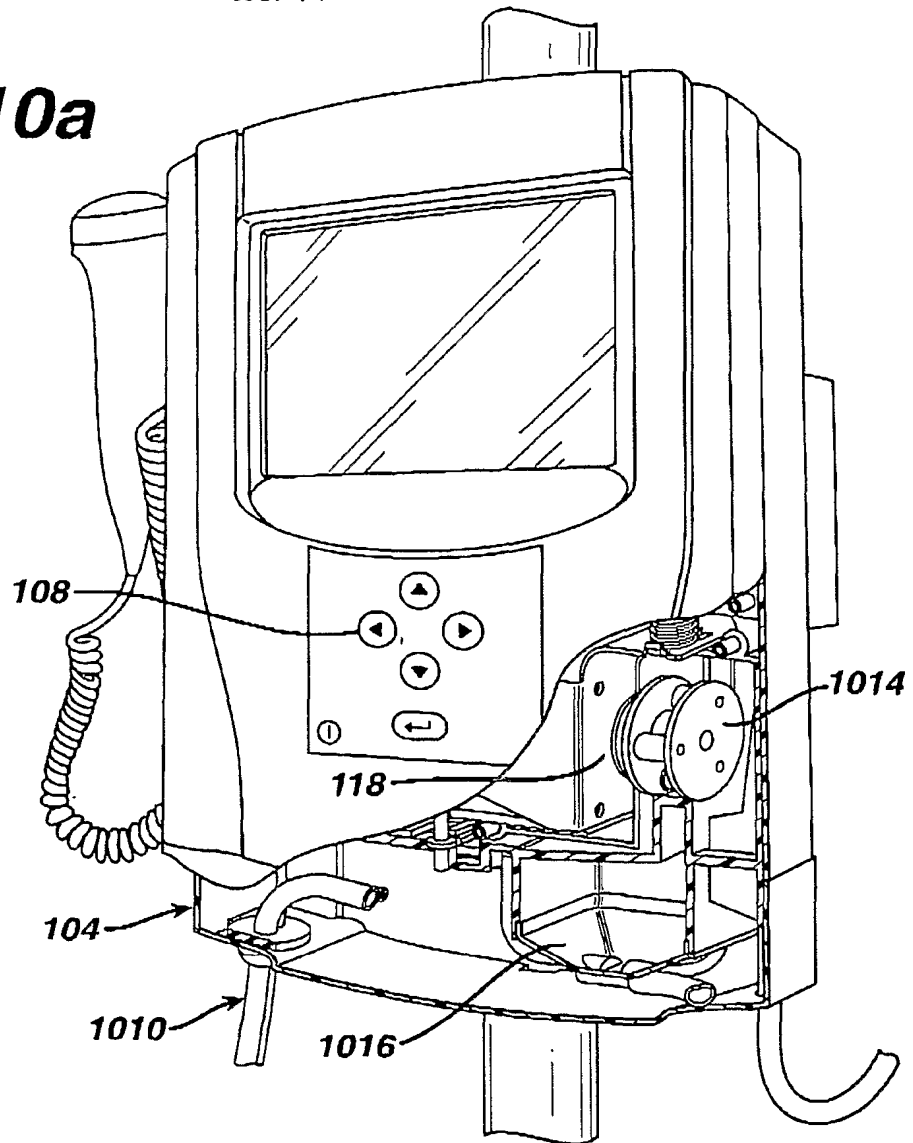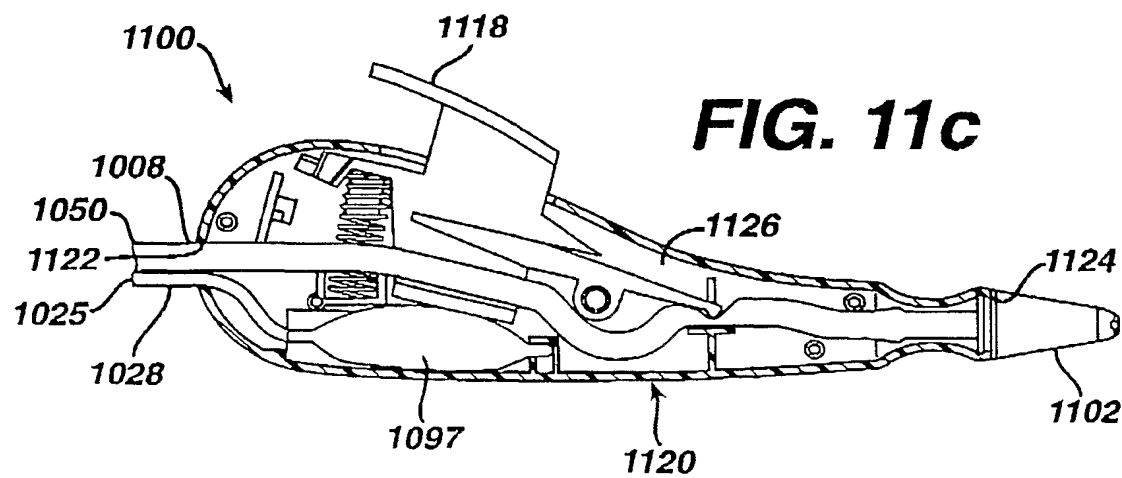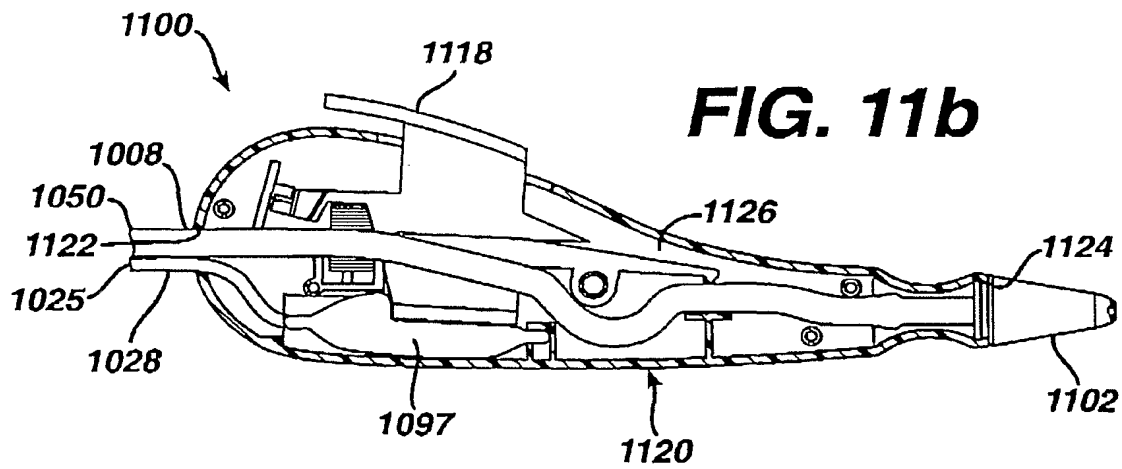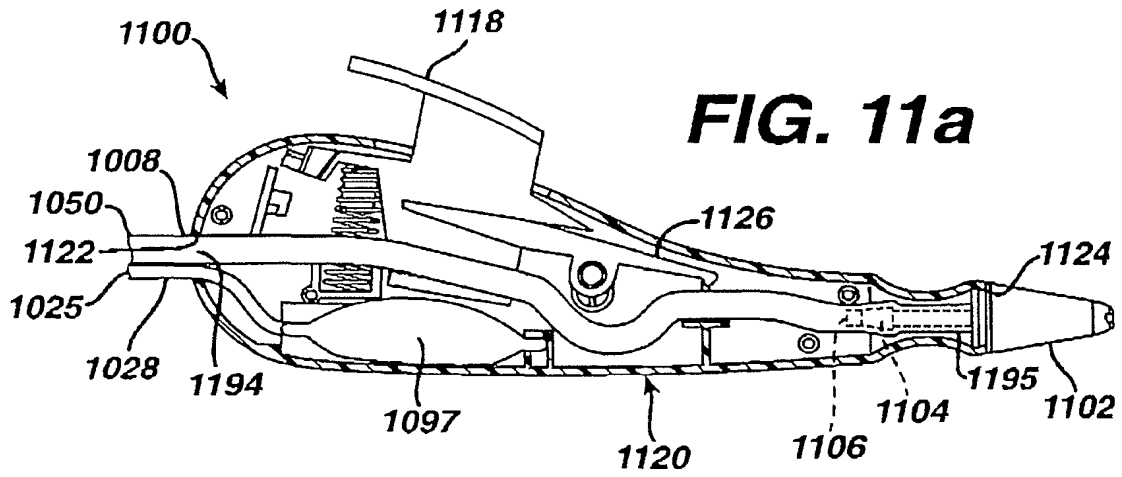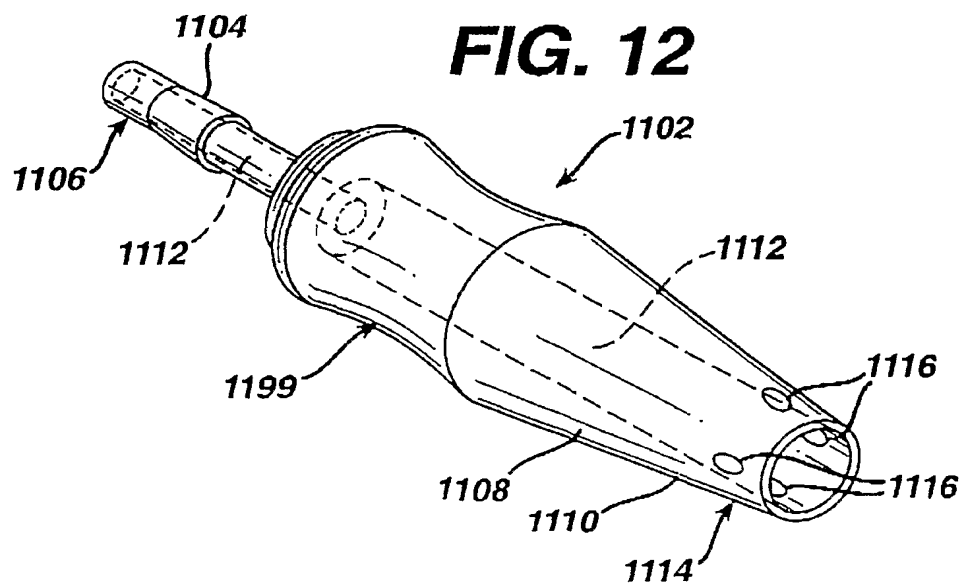

25/41

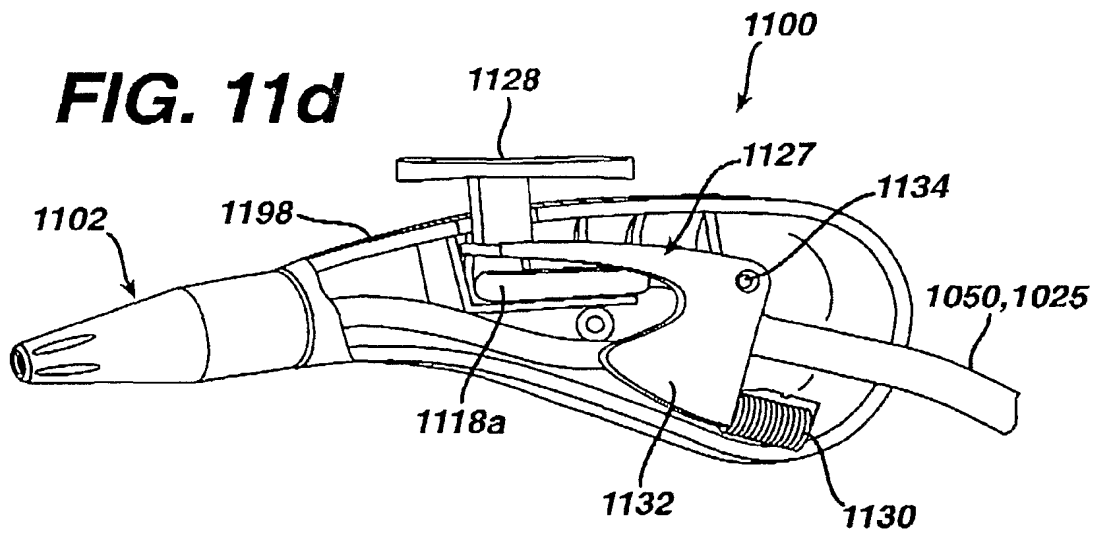

26/41



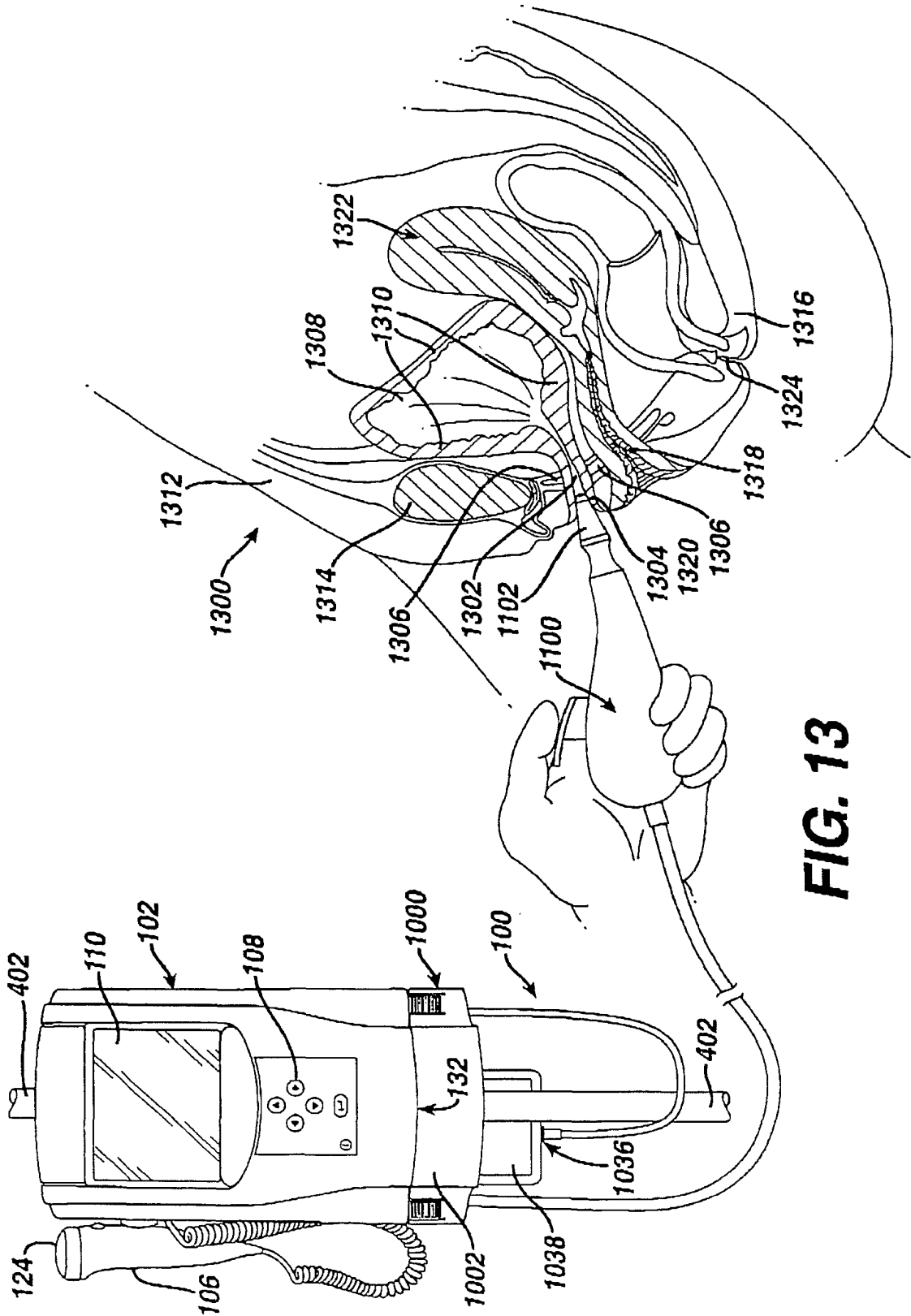



FIG. 13

27/41


FIG. 14

28/41

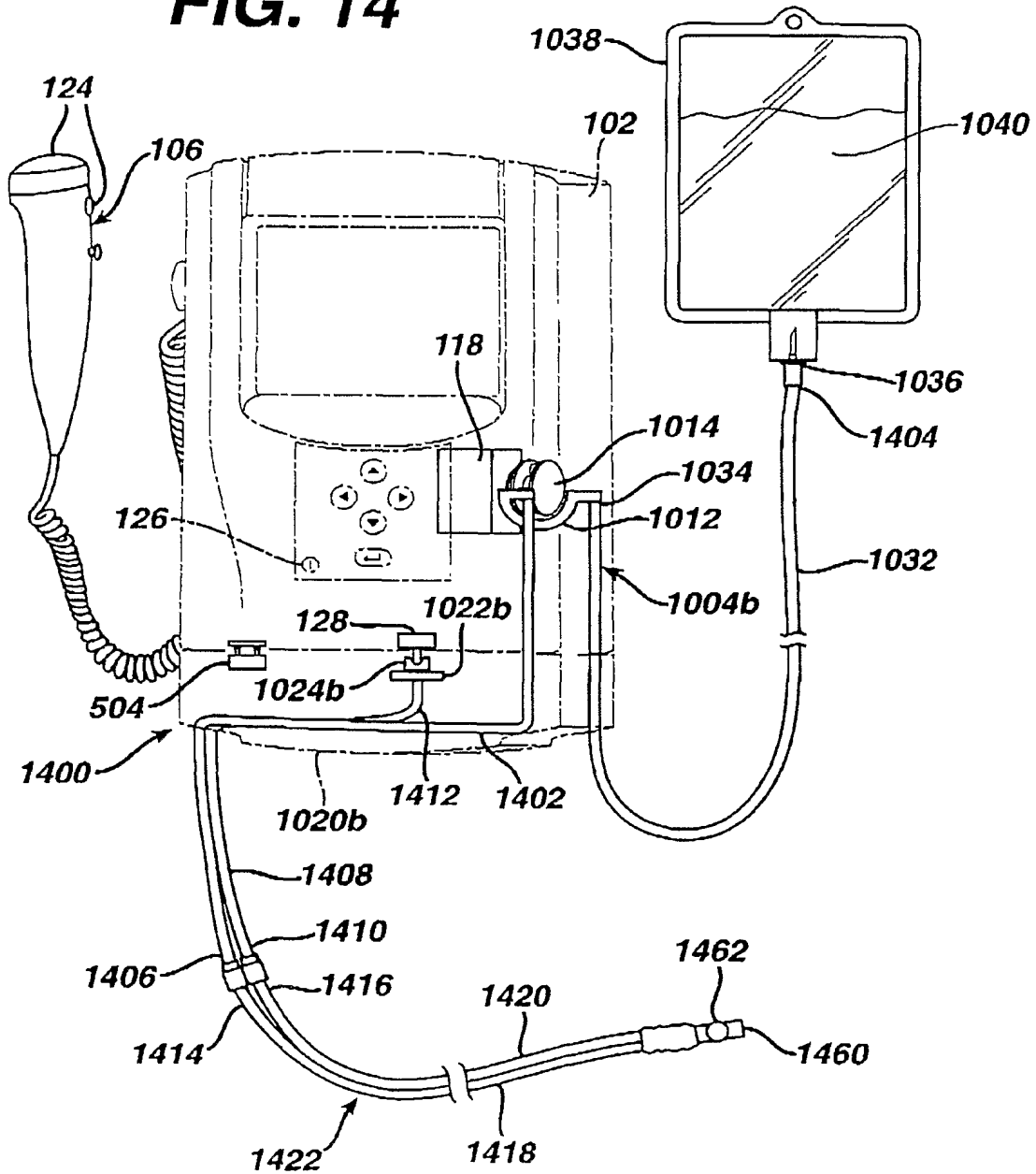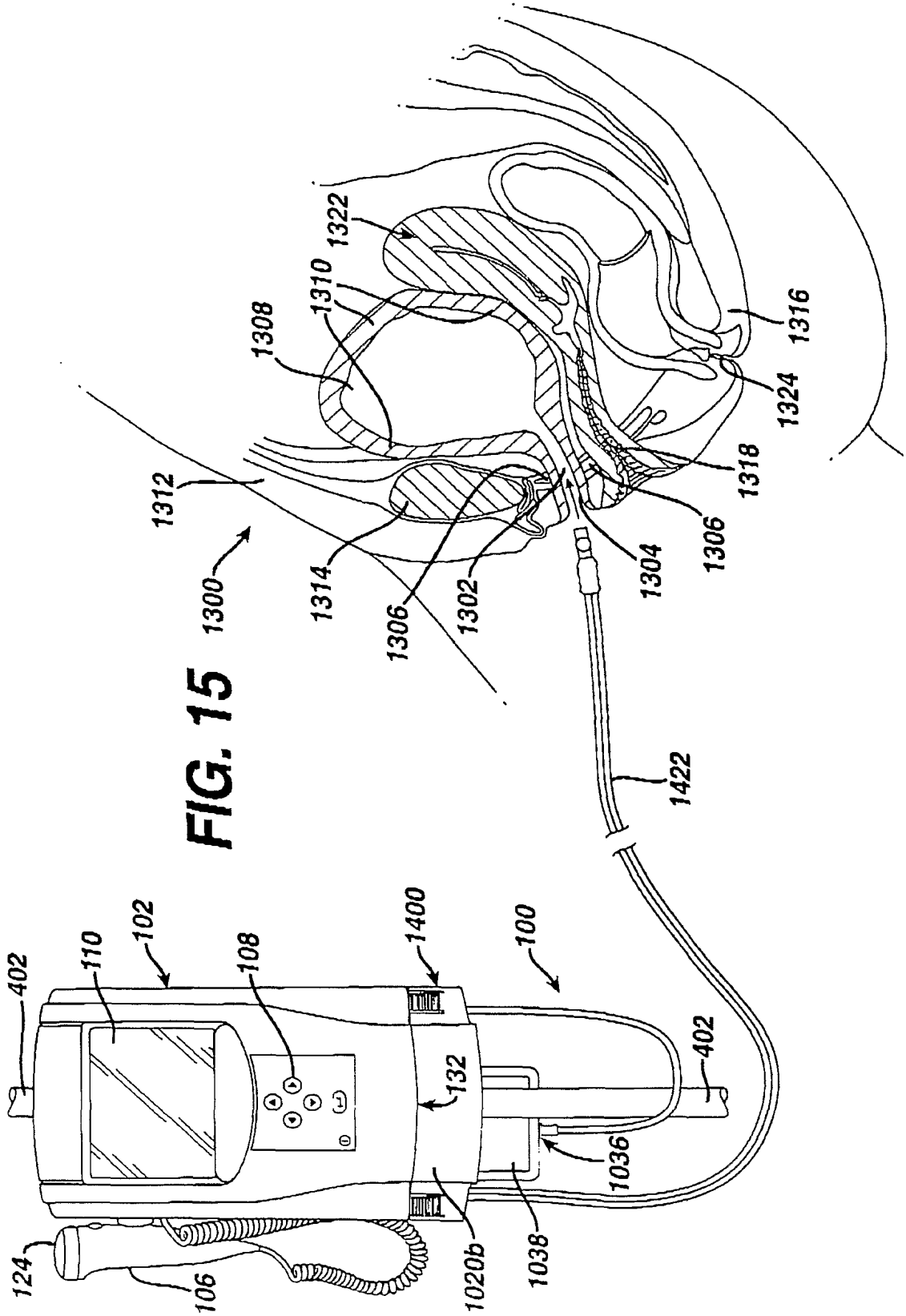


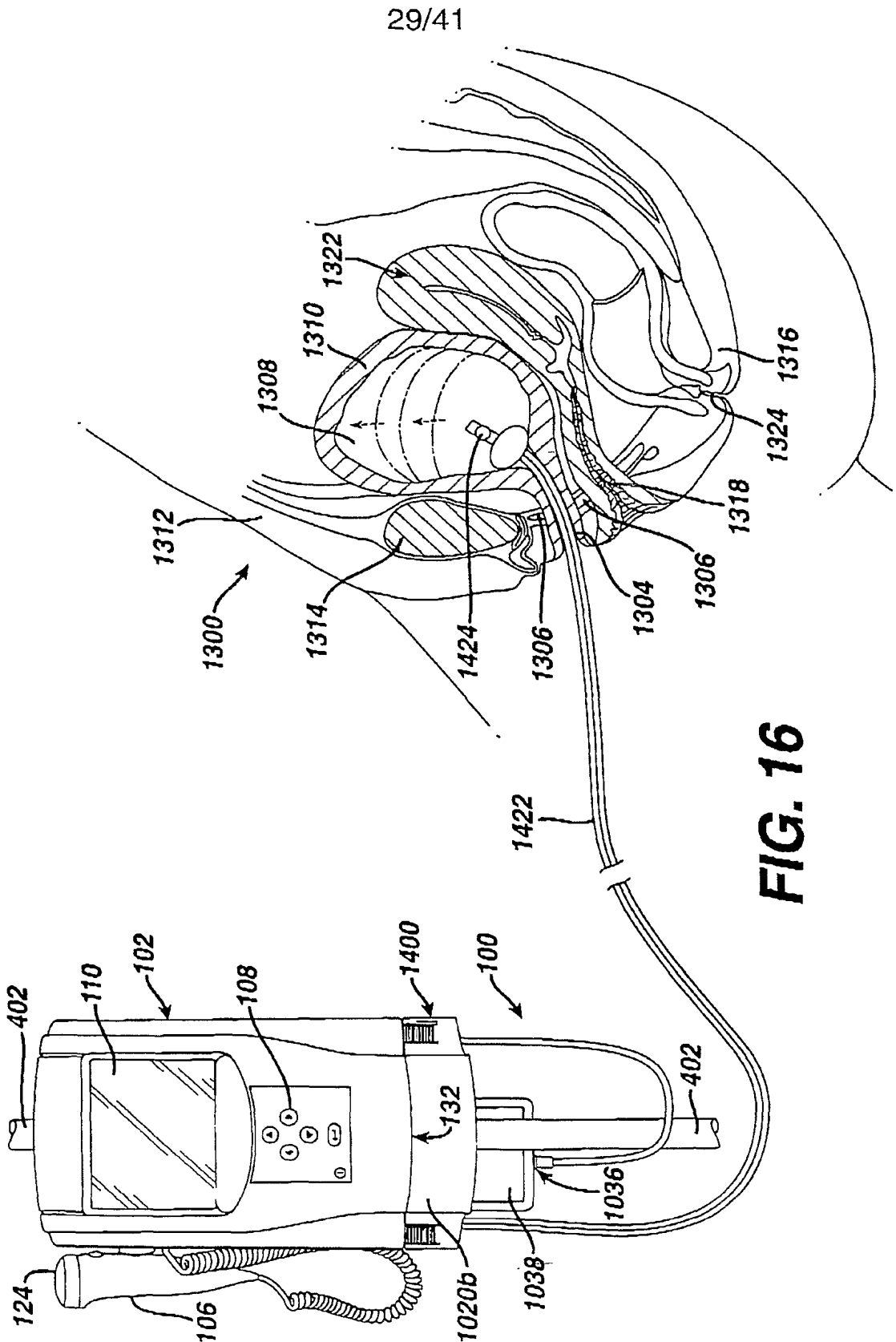


30/41

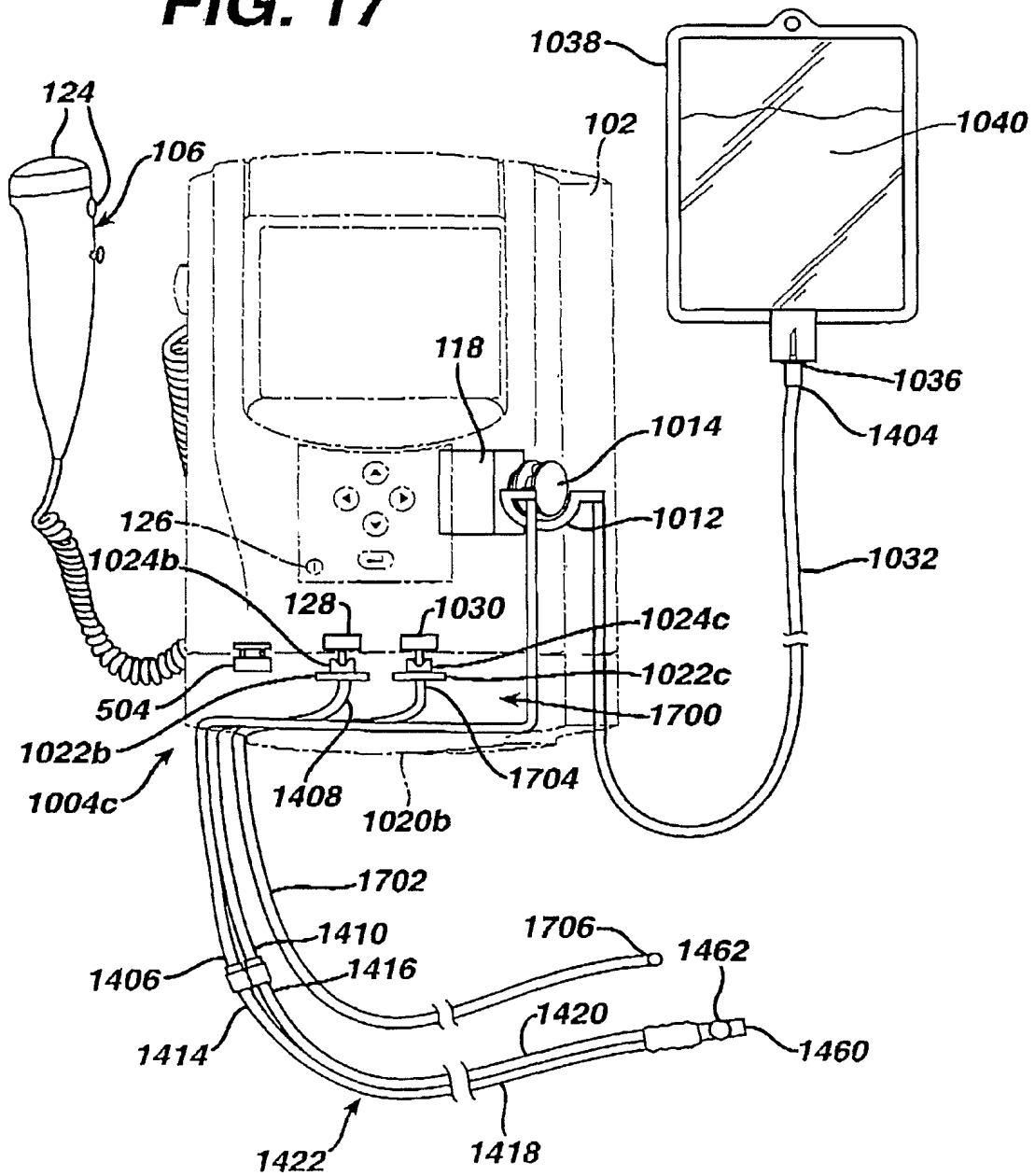FIG. 17


31/41

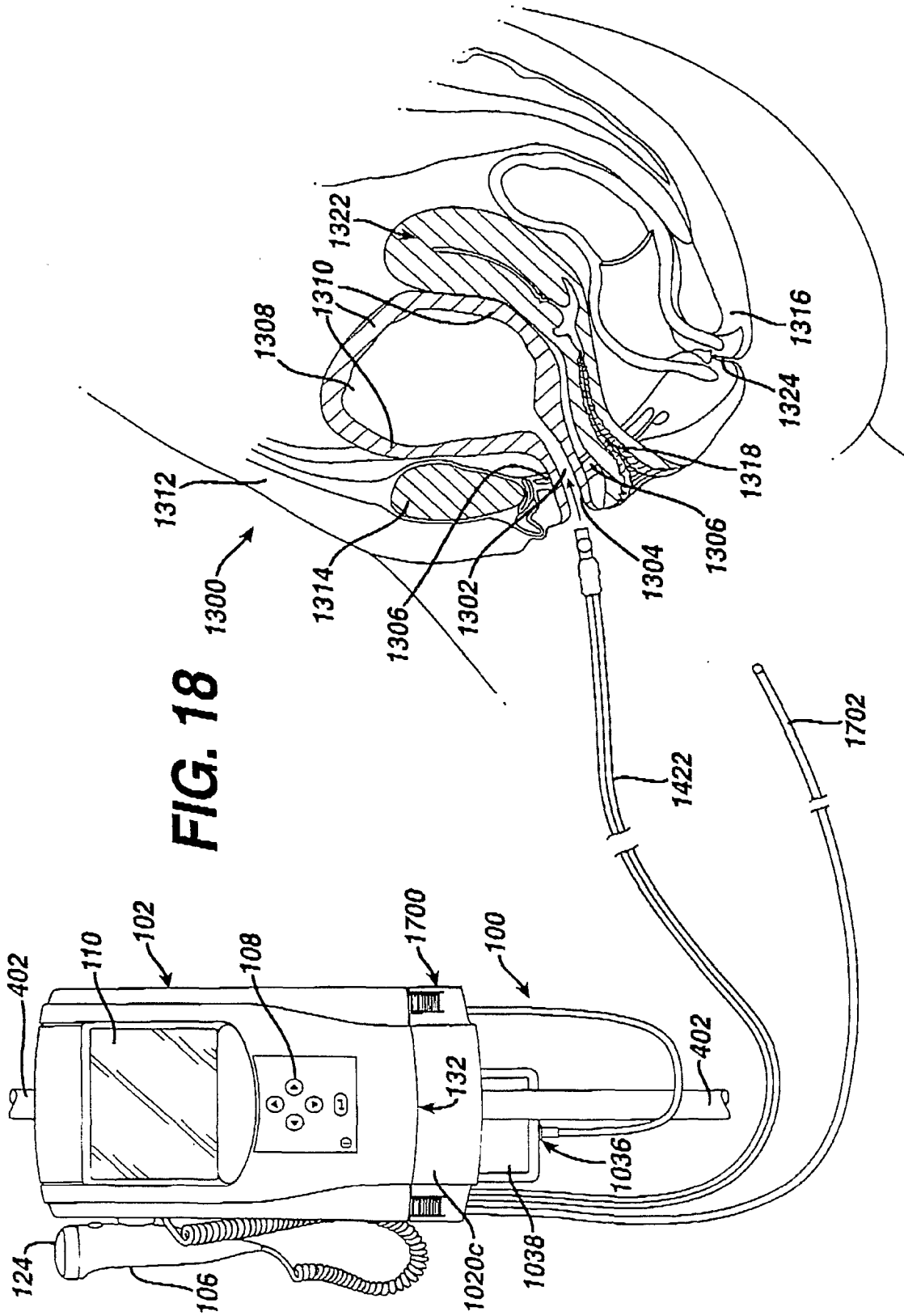


32/41

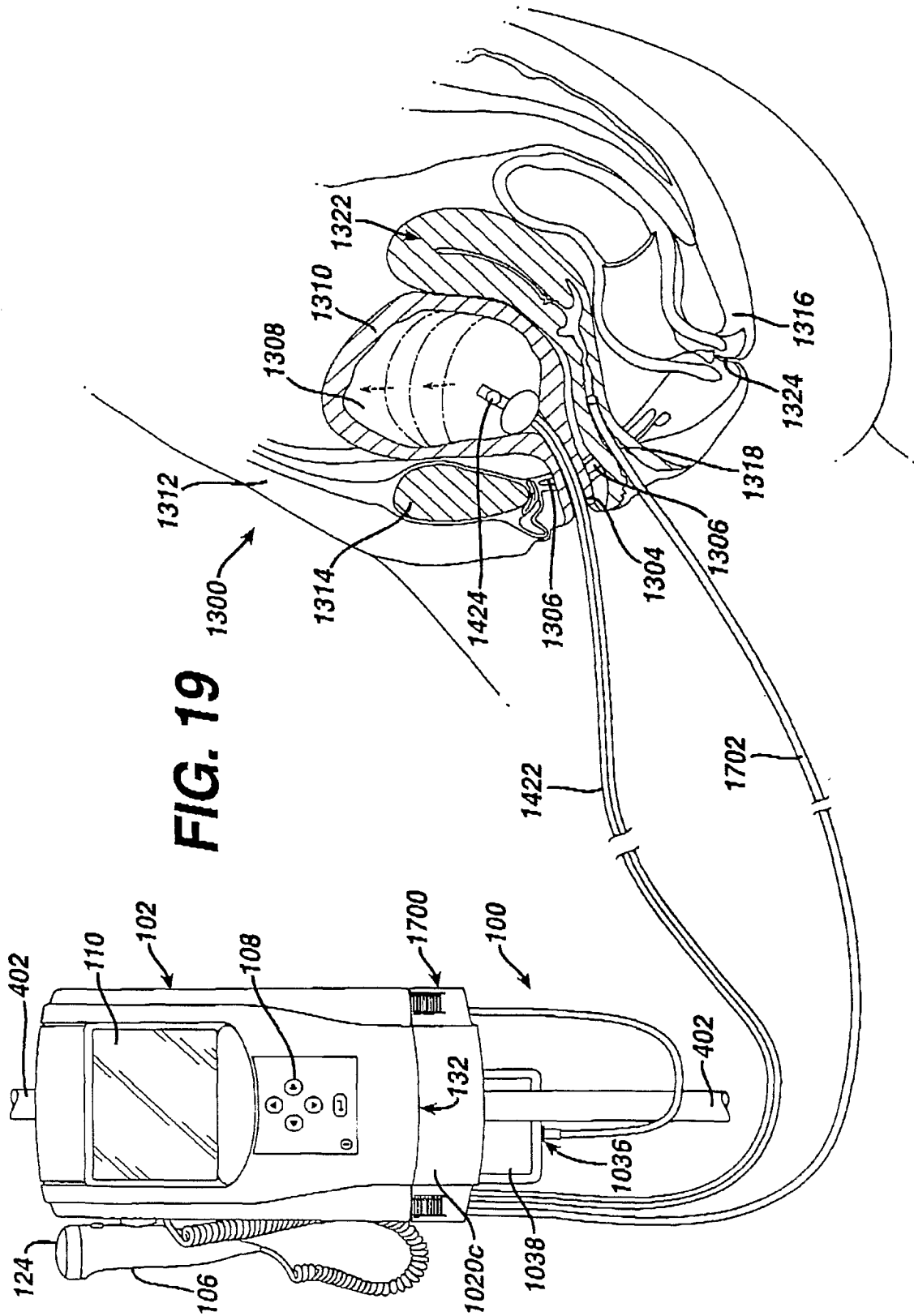


33/41

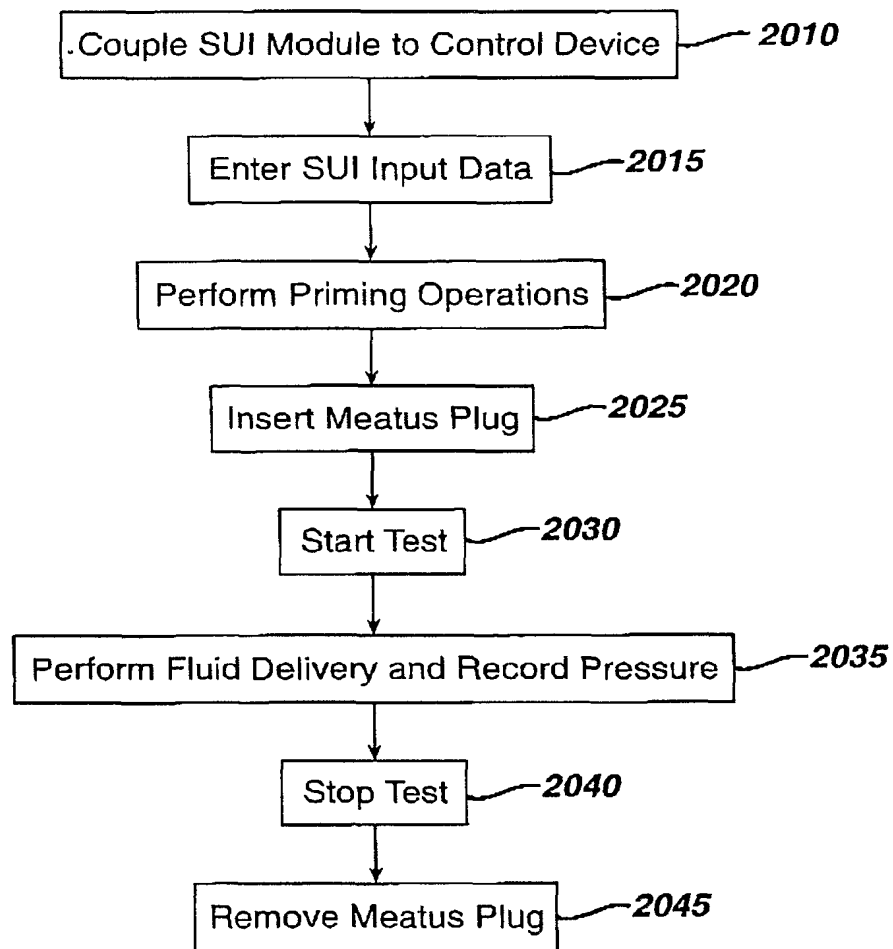FIG. 20

34/41

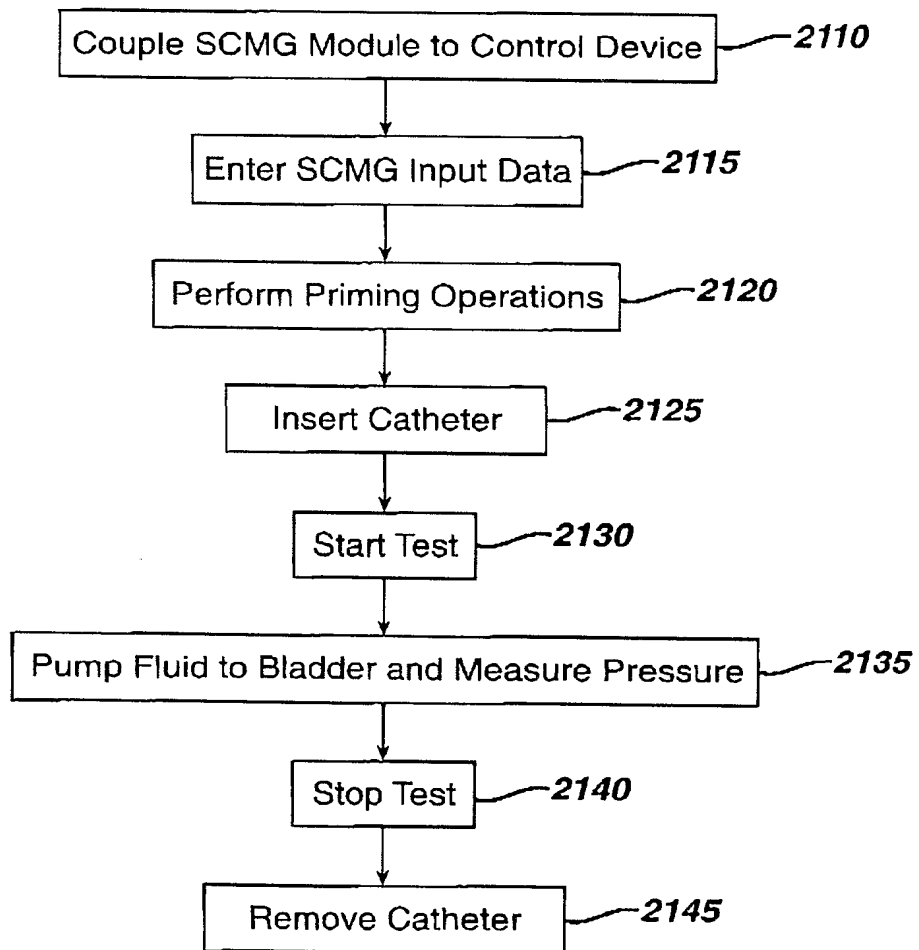FIG. 21

35/41

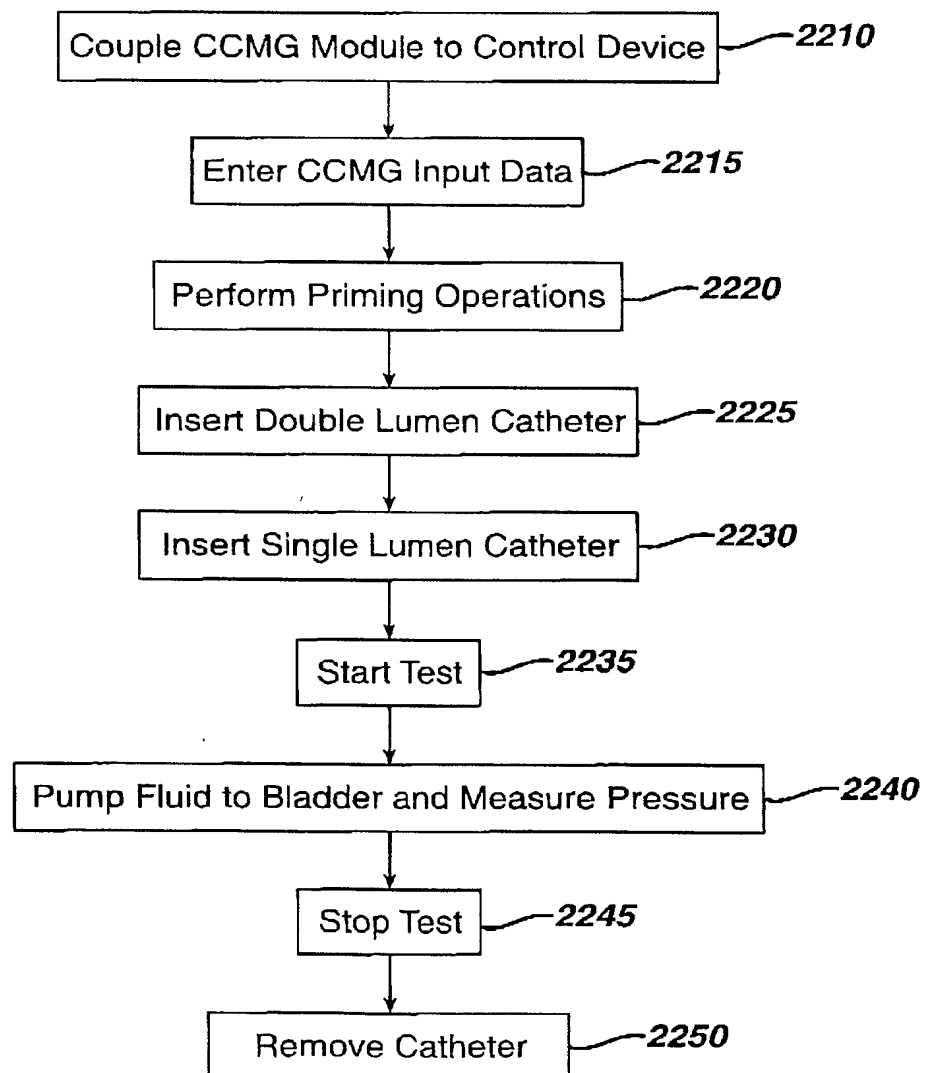FIG. 22

36/41

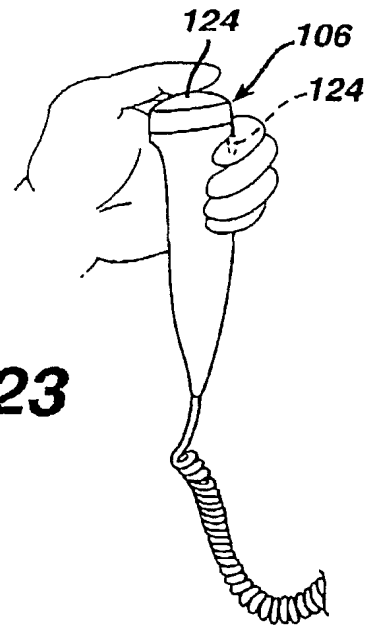


FIG. 23

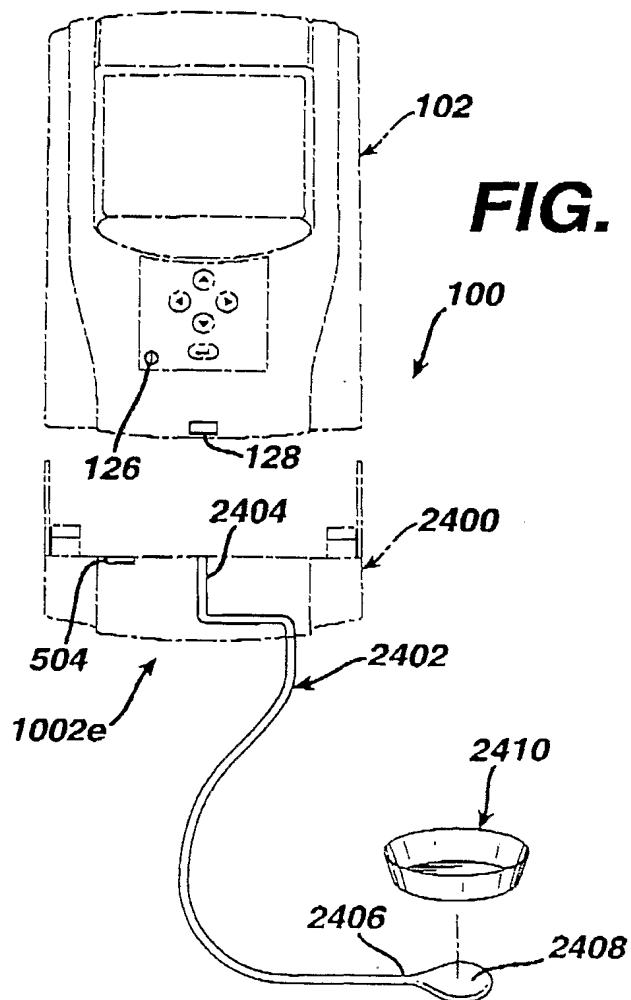
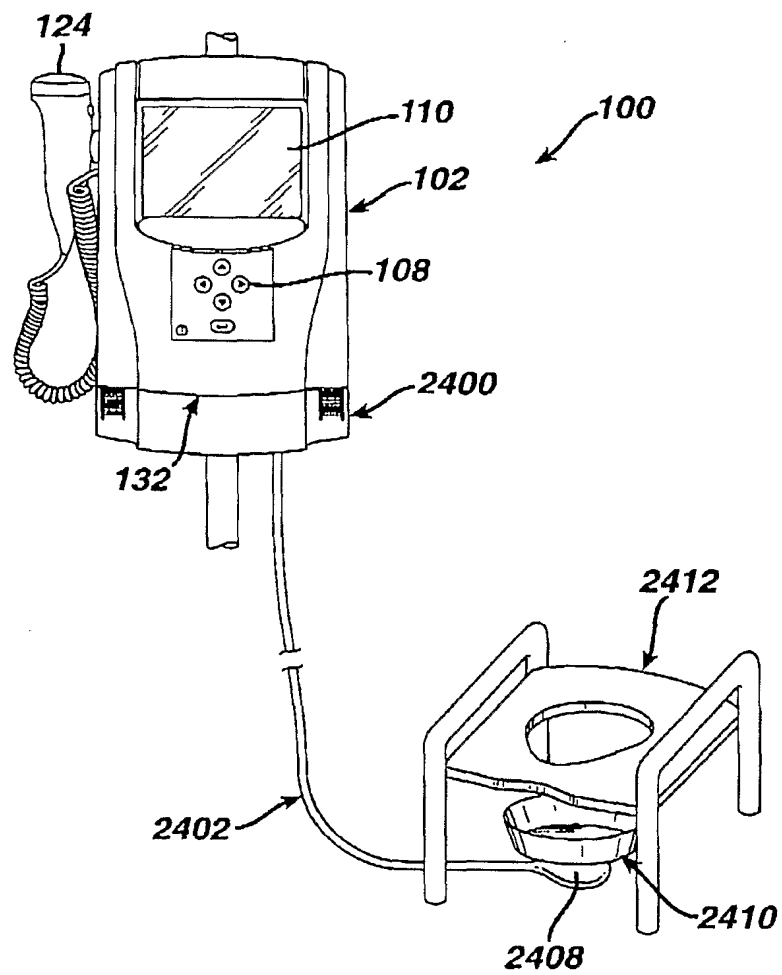
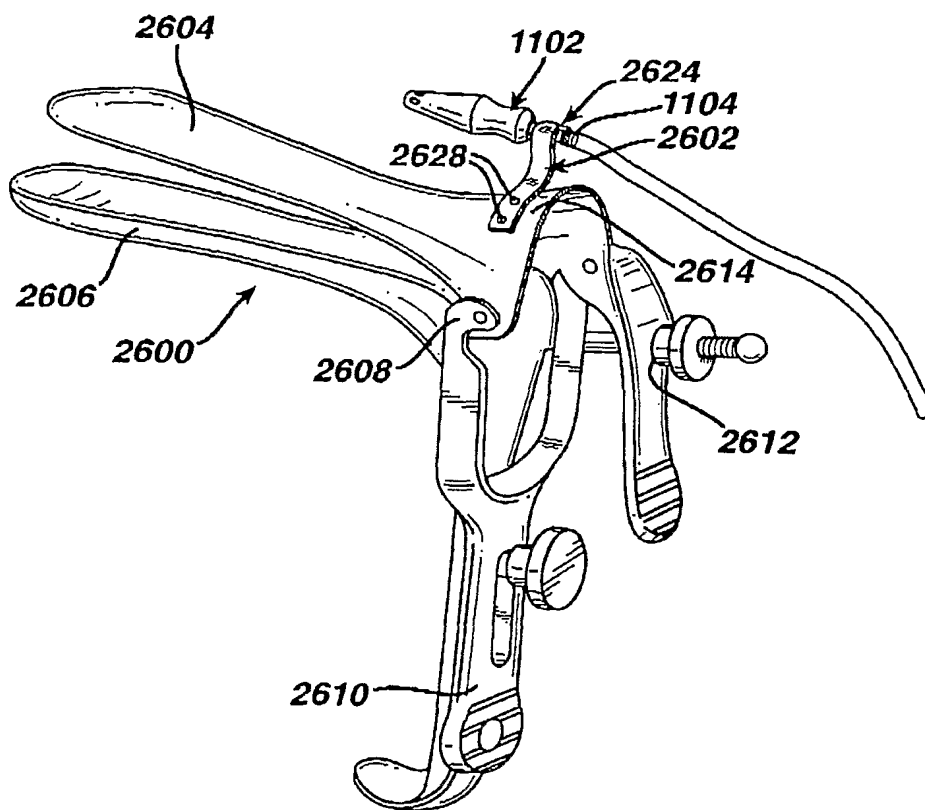


FIG. 24


37/41

FIG. 25


38/41

FIG. 26

39/41

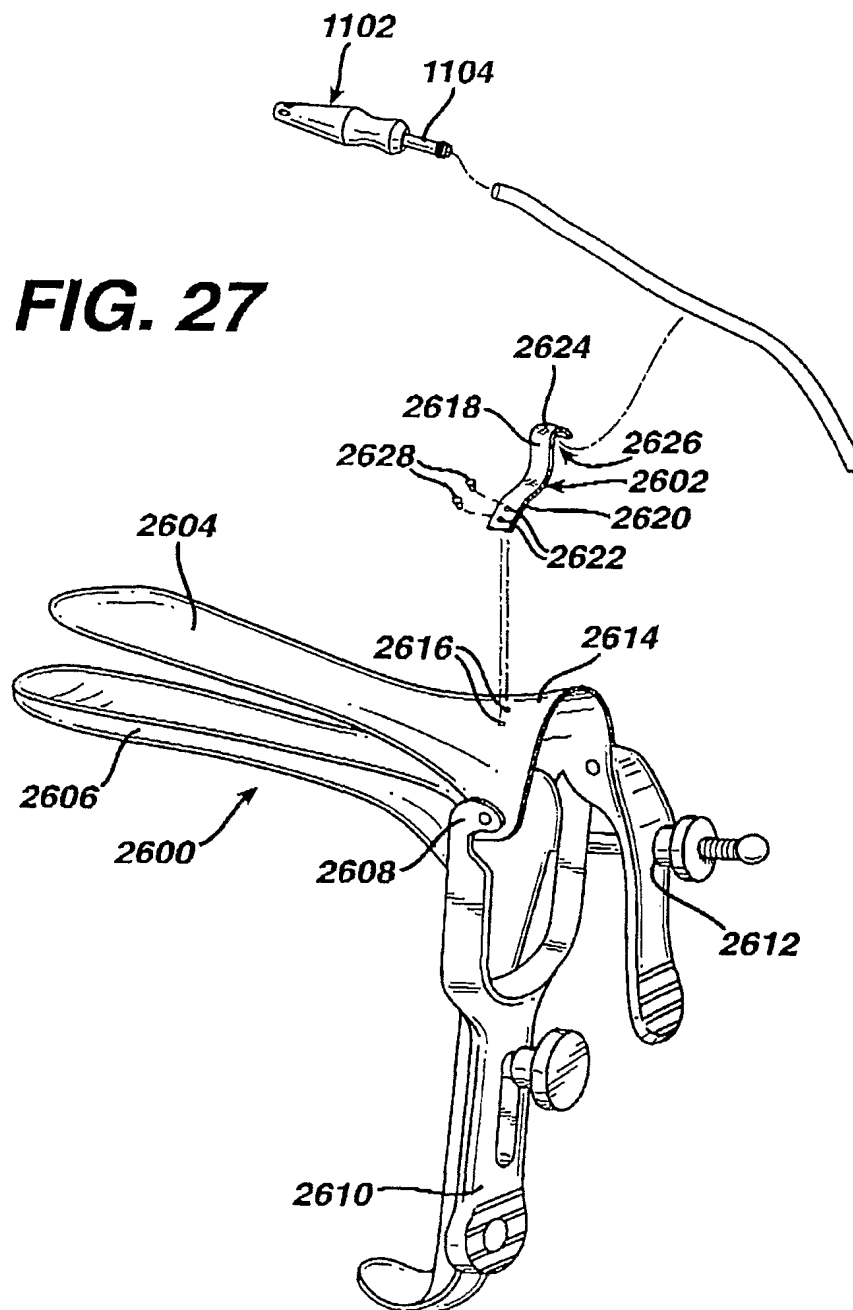

40/41

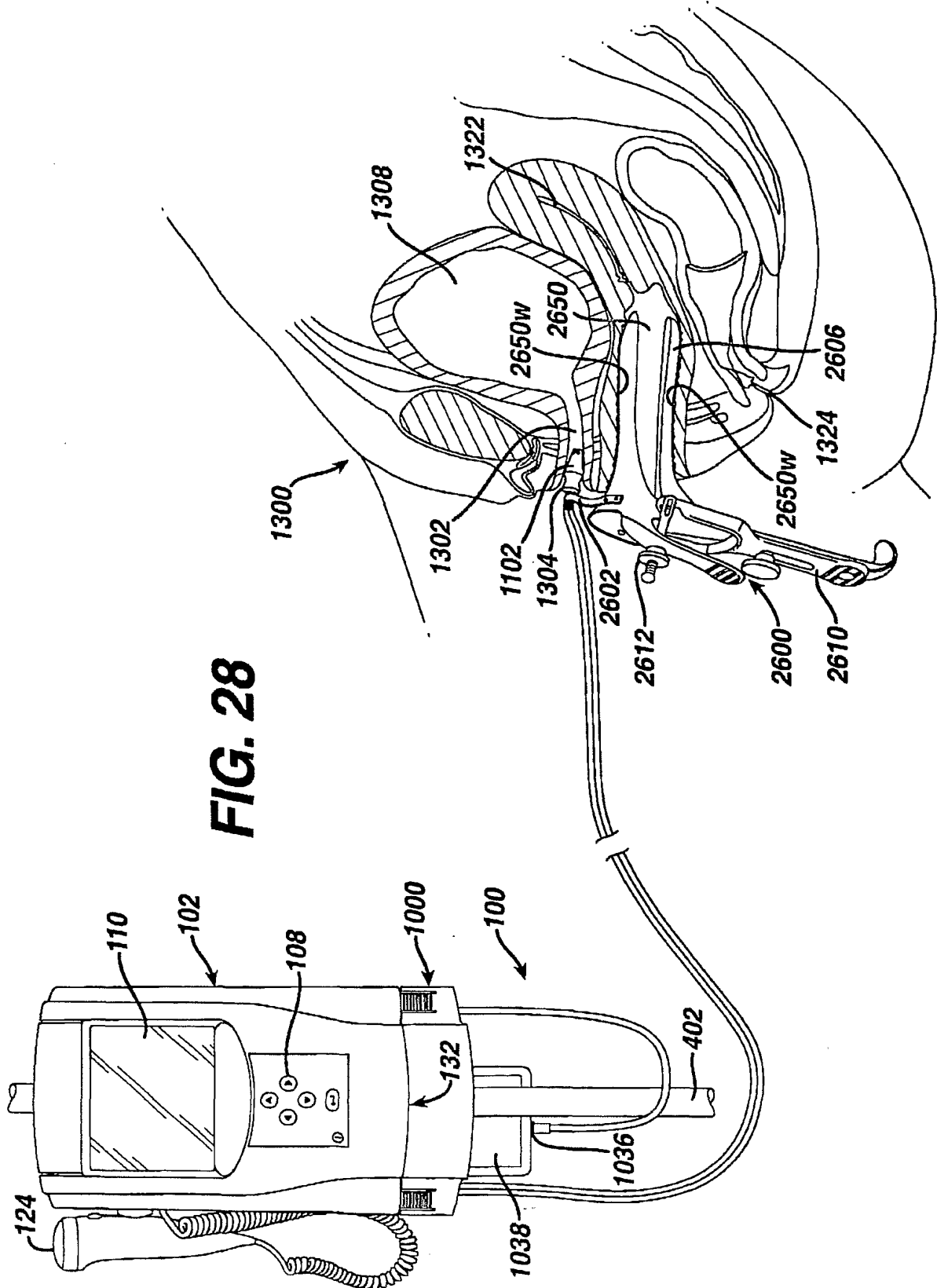

41/41

FIG. 29

